


Analysis of Nucleic acid base

Data No. LA816-0699

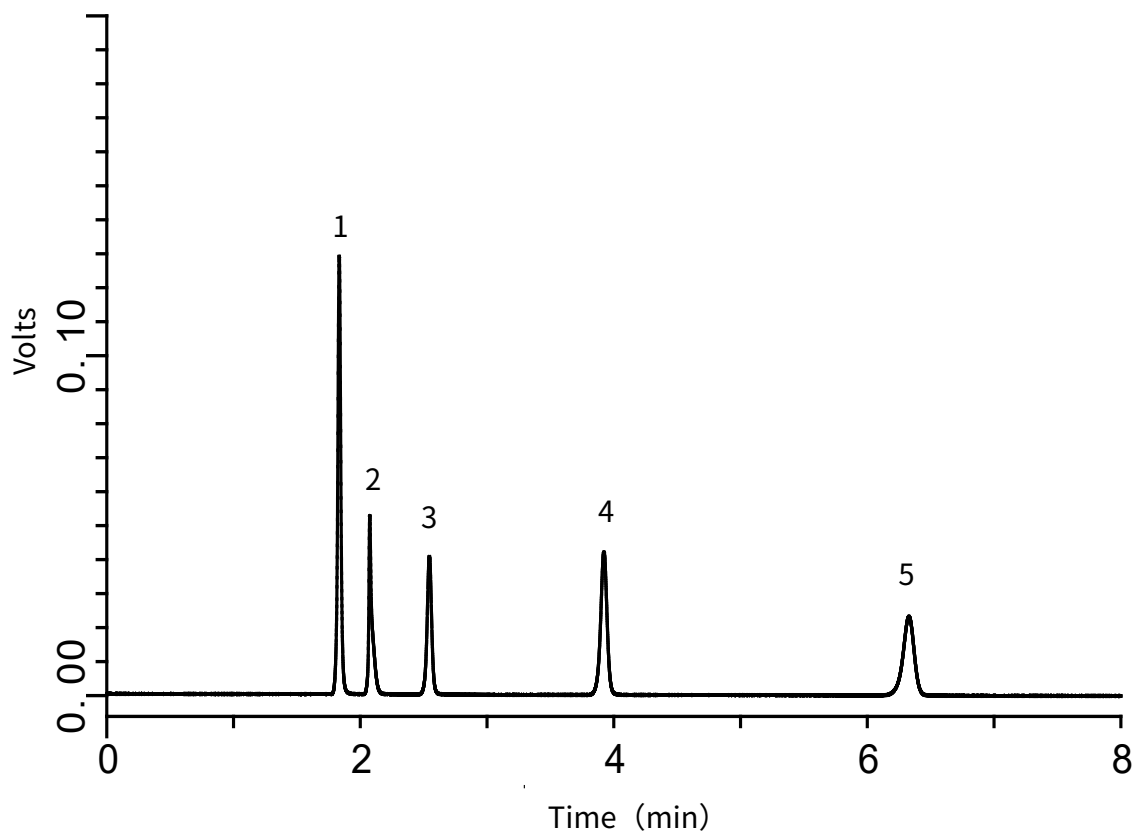

Conditions

Column : MonoClad C18-HS (250 x 3.0 mm I.D.)
Column Cat. No. : 5020-10804
Eluent : 0.1 M Na₂HPO₄ (pH 6.0, H₃PO₄)
Flow Rate : 1 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm
Injection Vol. : 0.5 µL

Analyte

1. Cytosine (0.2 mg/mL)
2. Uracil (0.1 mg/mL)
3. Cytidine (0.2 mg/mL)
4. Thymine (0.3 mg/mL)
5. Guanosine (0.2 mg/mL)