



InertSearch™ for LC

Inertsil® Applications

Analysis of Candesartan cilexetil

(Under the Condition of the Japanese Pharmacopoeia)

Data No. LB201-0919

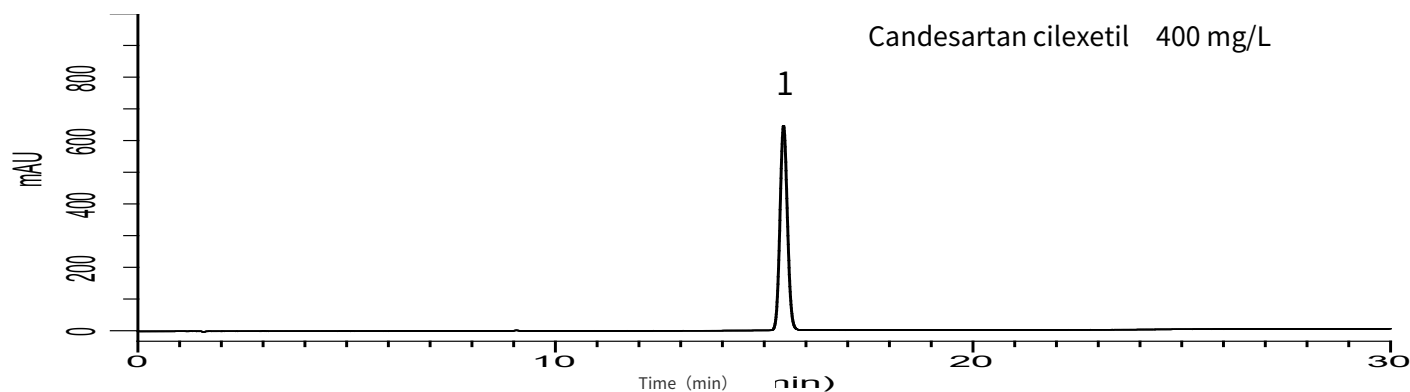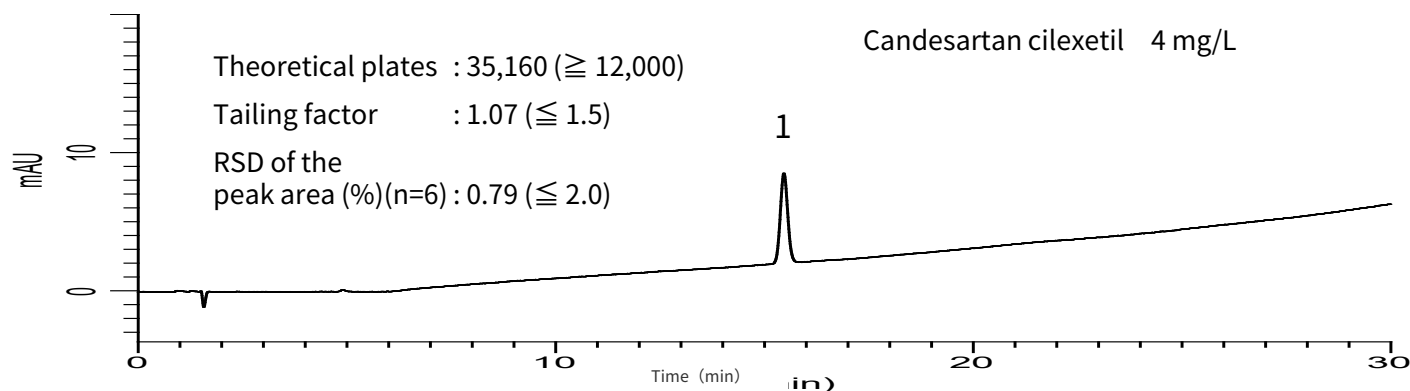

Conditions

System : GL-7400 HPLC system
Column : Inertsil ODS-3
(4 μ m, 150 x 4.0 mm I.D.)
Column Cat. No. : 5020-04635
Eluent : A) CH₃CN/H₂O/CH₃COOH = 57/43/1, v/v/v
B) CH₃CN/H₂O/CH₃COOH = 90/10/1, v/v/v
A/B = 100/0 - 30 min - 0/100, v/v
Flow Rate : 0.8 mL/min
Col. Temp. : 25 °C
Detection : UV 254 nm (GL-7452 PDA Detector)
Injection Vol. : 10 μ L
Sample : Standard

Analyte:

1. Candesartan cilexetil