


InertSearch™ for LC

Inertsil® Applications

Analysis of Vitamin E

Data No. LB168-0894

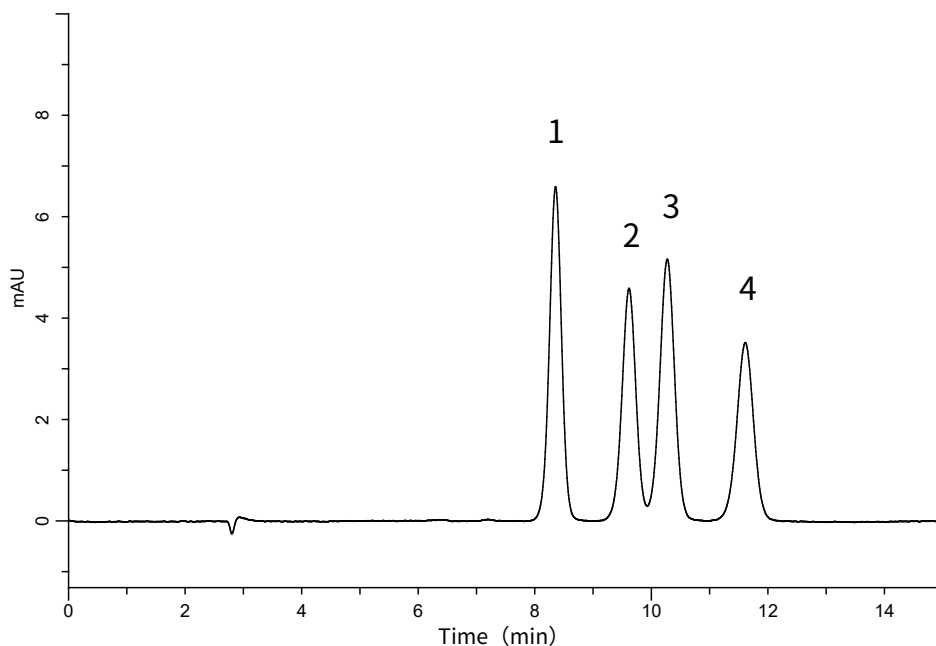

Conditions

System : GL-7400 HPLC system
Column : Inertsil C30 S-Select
(5 μ m, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-17037
Eluent : CH₃OH
Flow Rate : 1.0 mL/min
Col. Temp. : 15 °C
Detection : UV 295 nm (GL-7452 PDA Detector)
Injection Vol. : 1 μ L
Sample : Standard

Analyte:

1. δ -tocopherol	0.25 mg/mL
2. γ -tocopherol	0.25 mg/mL
3. β -tocopherol	0.25 mg/mL
4. α -tocopherol	0.25 mg/mL