


InertSearch™ for LC

Inertsil® Applications

Analysis of Resveratrol

Data No. LA522-0602

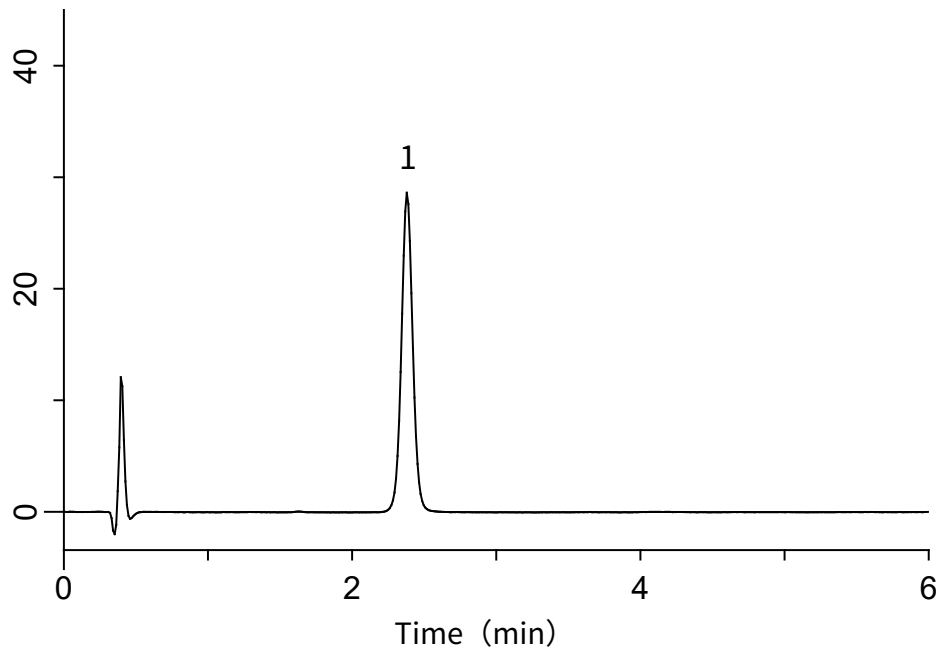

Conditions

System : GL7700 HPLC system
Column : InertSustain C18
(3 μ m, 50 x 3.0 mm I.D.)

Column Cat. No. : 5020-07422

Eluent : A) CH₃CN
B) H₂O
A/B = 20/80, v/v

Flow Rate : 1.0 mL/min

Col. Temp. : 40 °C

Detection : UV 310 nm (UV7750 UV Detector)

Injection Vol. : 20 μ L

Sample : Standard

Analyte:

1. Resveratrol 1.0 mg/L