


InertSearch™ for LC

Inertsil® Applications

Analysis of PAHs

Data No. LA766-0000

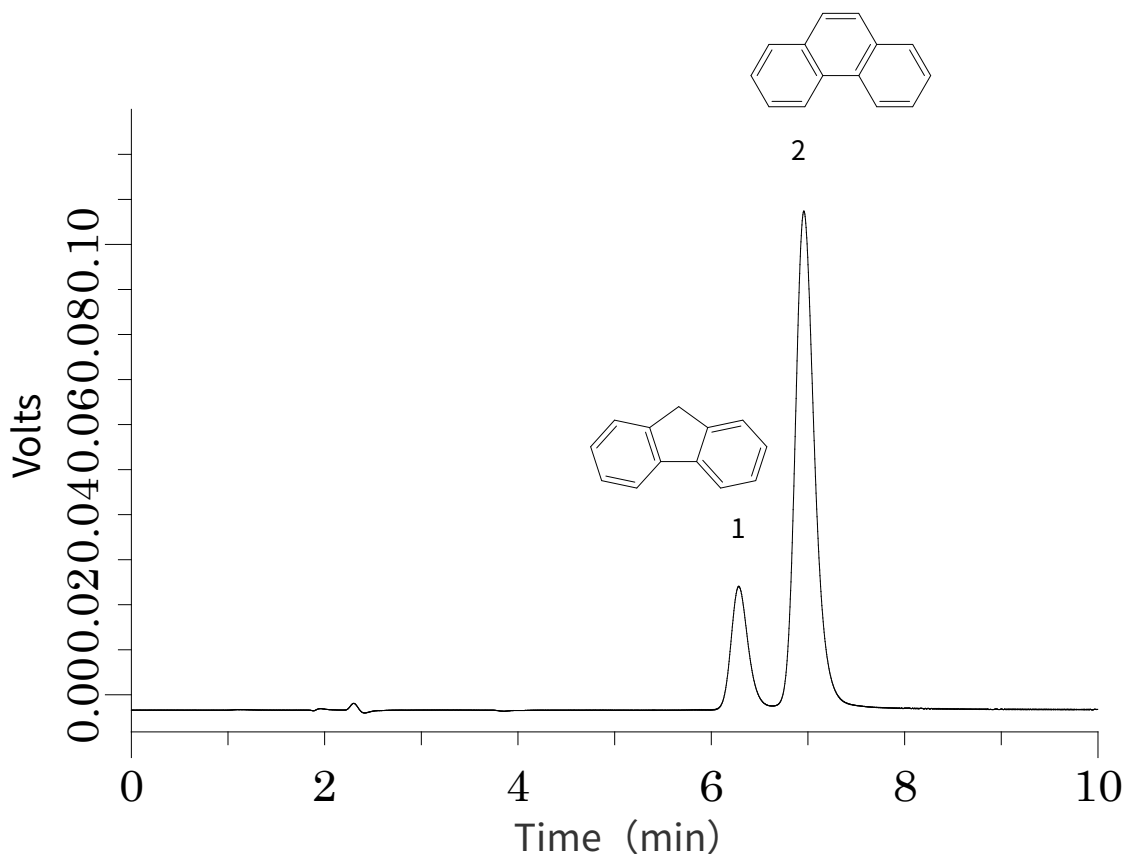

Conditions

Column : Inertsil Ph-3 (5 μ m, 150 x 4.6 mm I.D.)

Column Cat. No. : 5020-01920

Eluent : A) CH₃OH

B) H₂O

A/B = 70/30, v/v

Flow Rate : 1.0 mL/min

Col. Temp. : 30 °C

Detection : UV 254 nm

Analyte

1. Fluorene

2. Phenanthrene