


InertSearch™ for LC

Inertsil® Applications

Analysis of Recerpine (Inertsil ODS-4)

Data No. LA976-0811

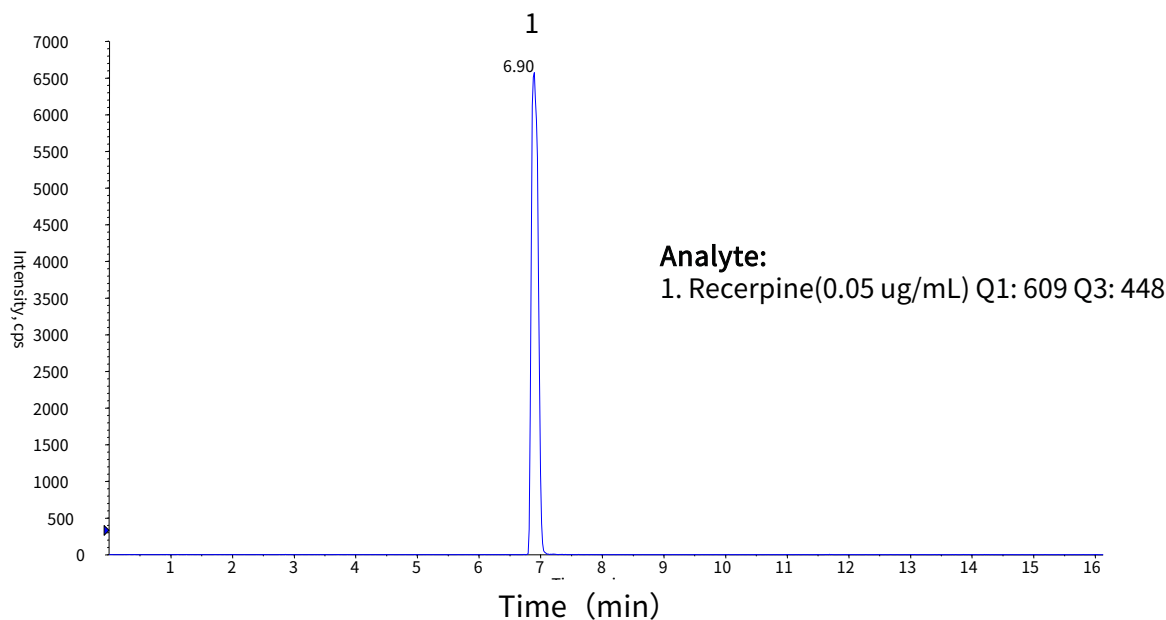

Conditions

System	: LC800 HPLC system 4000 Q TRAP®
Column	: Inertsil ODS-4 HP (3µm, 50 x 2.1 mm I.D.)
Column Cat. No.	: 5020-14062
Eluent	: A) 0.02 % TFA in CH ₃ CN : B) 0.02 % TFA in H ₂ O A/B = 25/75 –3 min- 25/75 –2 min- 90/10 -7 min- 90/10, v/v
Flow rate	: 0.2 mL/min
Col. Temp.	: 40 °C
Detection	: LC/MS/MS (4000Q TRAP® : ESI, Positive, MRM) CUR CAD IS TEM GS1 GS2 20 6 4500 600 70 80
Injection Vol.	: 10 µL
Sample	: Standard