


## Analysis of Tocopherol acetate

Data No. LA802-0696

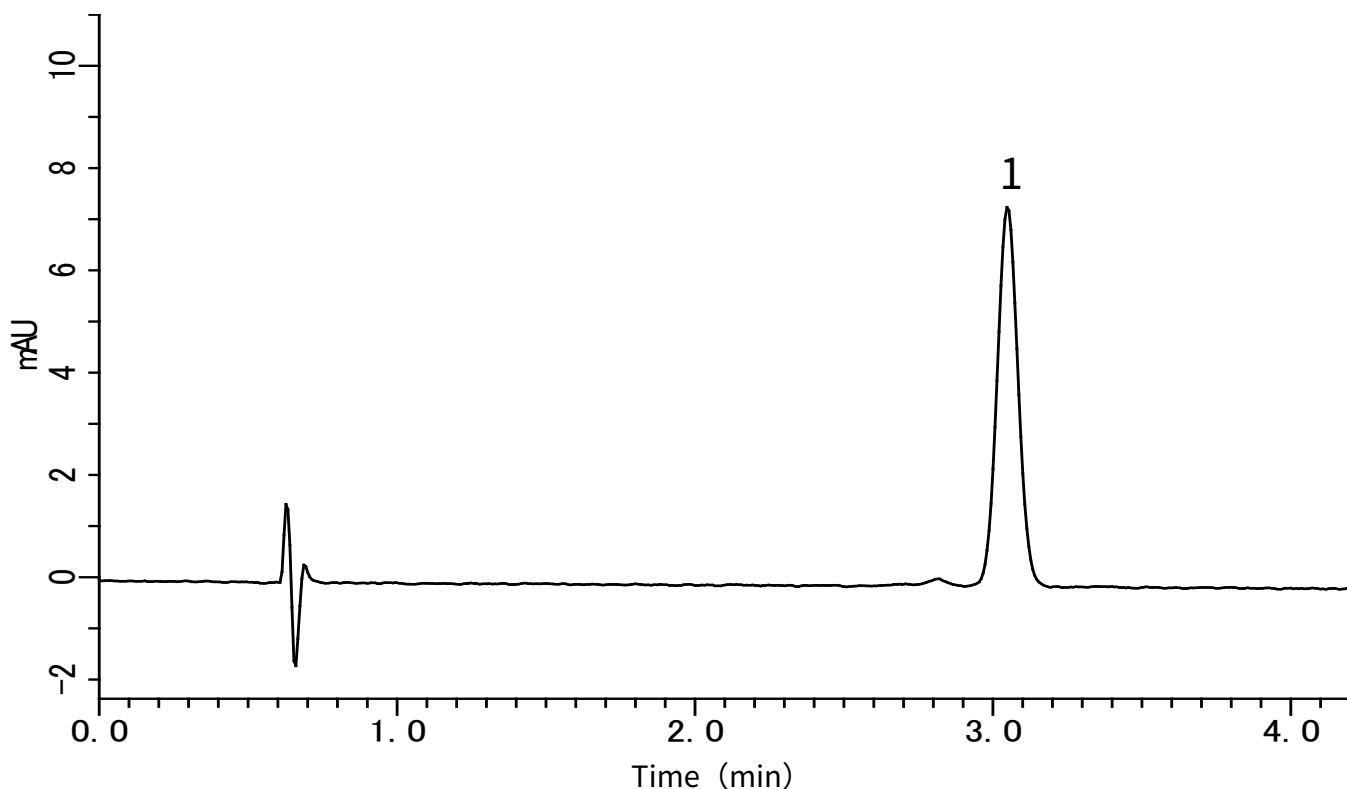

### Conditions

**System** : GL-7400 HPLC system  
**Column** : Inertsil ODS-4 (3μm, 100 x 3.0 mm I.D.)  
**Column Cat. No.** : 5020-04024  
**Eluent** : A) CH<sub>3</sub>CN  
          : B) CH<sub>3</sub>OH  
          : A/B = 80/20, v/v  
**Flow rate** : 0.8 mL/min  
**Col. Temp.** : 40 °C  
**Detection** : UV 285 nm (GL-7452A PDA Detector)  
**Injection Vol.** : 5 μL  
**Sample** : Tocopherol acetate

### **Analyte:**

1. Tocopherol acetate 50 mg/L