


InertSearch™ for LC

Inertsil® Applications

Analysis of Tocopherol

Data No. LB227-0785

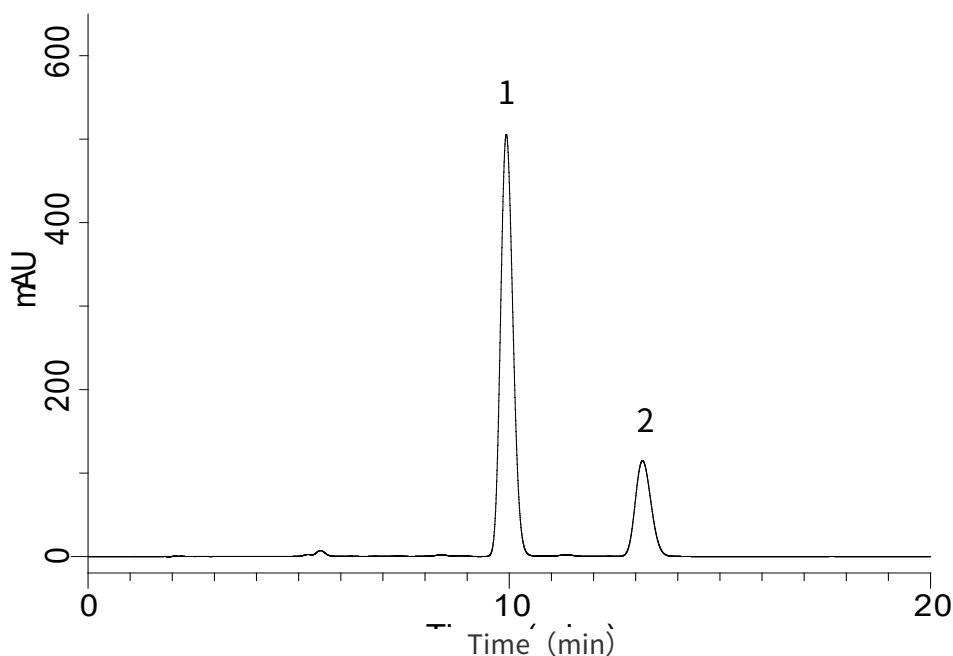

Conditions

System : GL7700 HPLC system
Column : Inertsil ODS-2
(5 μ m, 150 x 4.0 mm I.D.)
Column Cat. No. : 5020-01123
Eluent : A) CH₃OH
B) H₂O
A/B = 98/2, v/v
Flow Rate : 0.73 mL/min
Col. Temp. : 30 °C
Detection : UV 292 nm (UV7750 UV Detector)
Injection Vol. : 20 μ L
Sample : Standard

Analyte:

1. <i>d</i> - α -Tocopherol	1.0 mg/L
2. Tocopherol acetate	1.0 mg/L