


Analysis of Peptides and proteins (WP300 C4)

Data No. LA775-0000

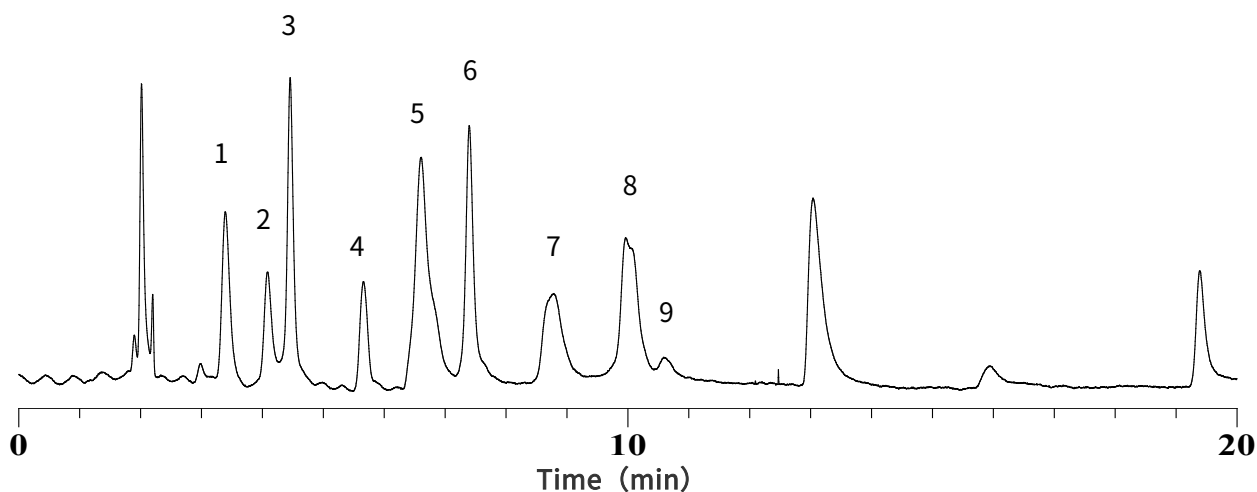

Conditions

Column : Inertsil WP300 C4 (5µm, 150 x 4.6 mm I.D.)

Column Cat. No. : 5020-05895

Eluent : A) 0.2 % HCOOH in (CH₃CN/H₂O = 90/10, v/v)

B) 0.2 % HCOOH

A/B = 20/80 – 20 min – 80/20, v/v

Flow Rate : 1.0 mL/min

Col. Temp. : 30 °C

Detection : UV 280 nm

Analyte

- | | |
|------------------------------|--------------|
| 1. Neurotensin | (FW 1673) |
| 2. Leucin Enkephalin | (FW 556) |
| 3. Cytochrome C | (FW 12,000) |
| 4. Insulin | (FW 6,000) |
| 5. BSA | (FW 66,000) |
| 6. Myoglobin | (FW 17,000) |
| 7. Creatine amidinohydrolase | (FW 43,000) |
| 8. Ovalbumin | (FW 45,000) |
| 9. Creatinine amidohydrolase | (FW 170,000) |