


InertSearch™ for LC

Inertsil® Applications

Analysis of Nucleotides

Data No. LA690-0000

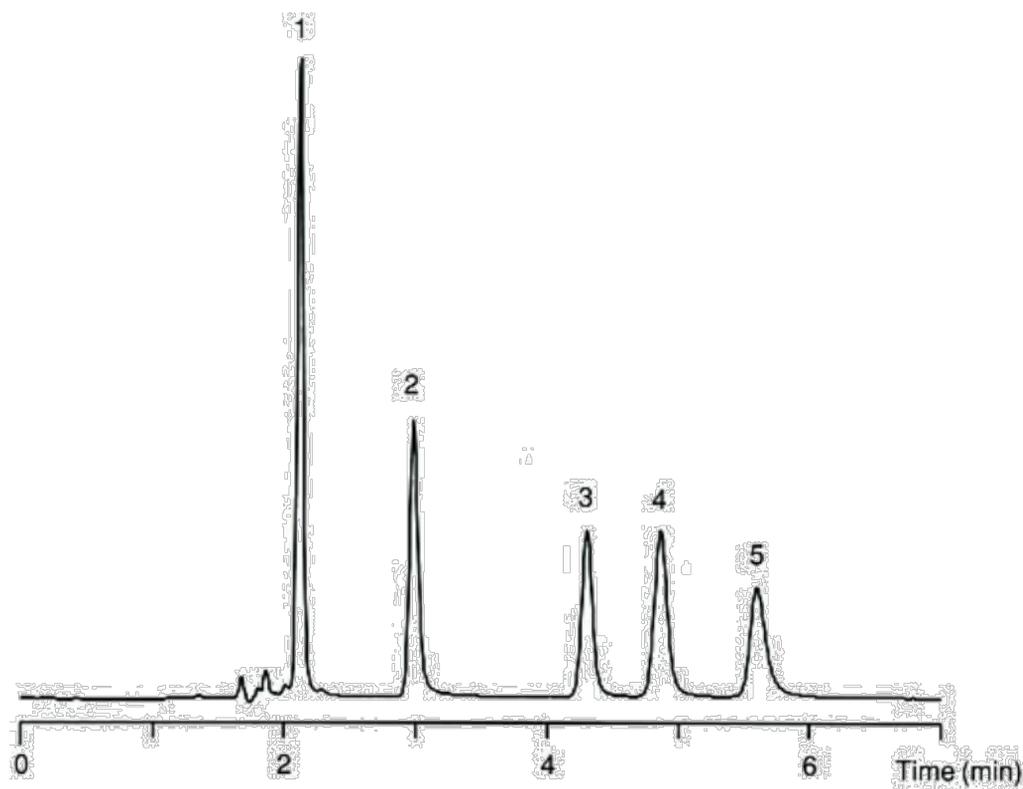

Conditions

Column : Inertsil NH2 (5 μ m, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-05545
Eluent : 25 mM Phosphate bufer (pH 2)
Flow Rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm

Analyte

1. 5' -AMP
2. 5' -GMP
3. 5' -UMP
4. 5' -IMP
5. 5' -TMP