


# InertSearch™ for LC

Inertsil® Applications

## Analysis of Surfactant

Data No. LA699-0734

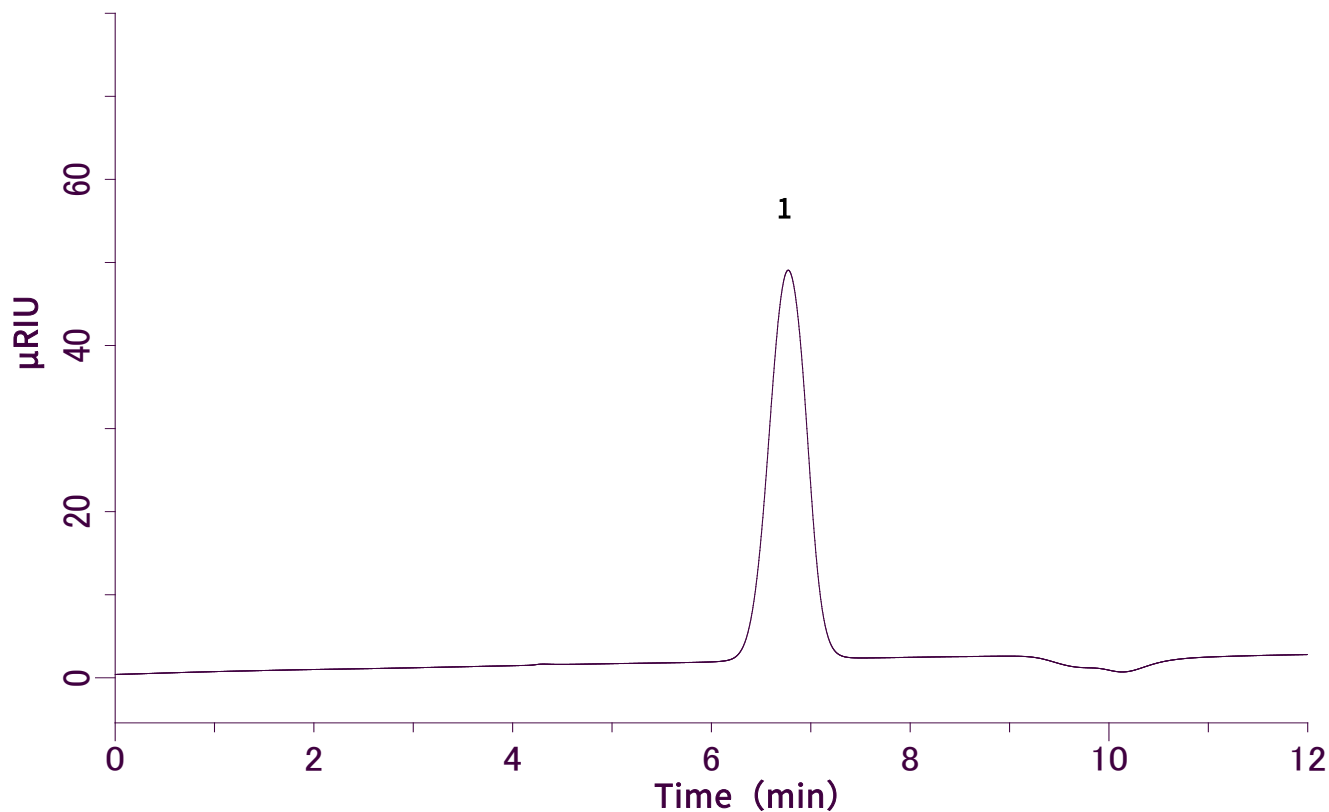

### Conditions

**System** : GL7400 HPLC system  
**Column** : Inertsil Diol (5μm, 250 x 7.6 mmI.D.)  
**Column Cat. No.** : 5020-05666  
**Eluent** : THF  
**Flow Rate** : 1.0 mL/min  
**Col. Temp.** : 40°C  
**Detection** : RI (GL-7454 RI Detector) (Cell Temp. 35°C)  
**Injection Vol.** : 50μL

### Analyte

1. Poly (oxyethylene) octylphenyl ether (Triton X-114)