


InertSearch for GC

InertCap Applications

Pesticides

Data No. GA321-1030

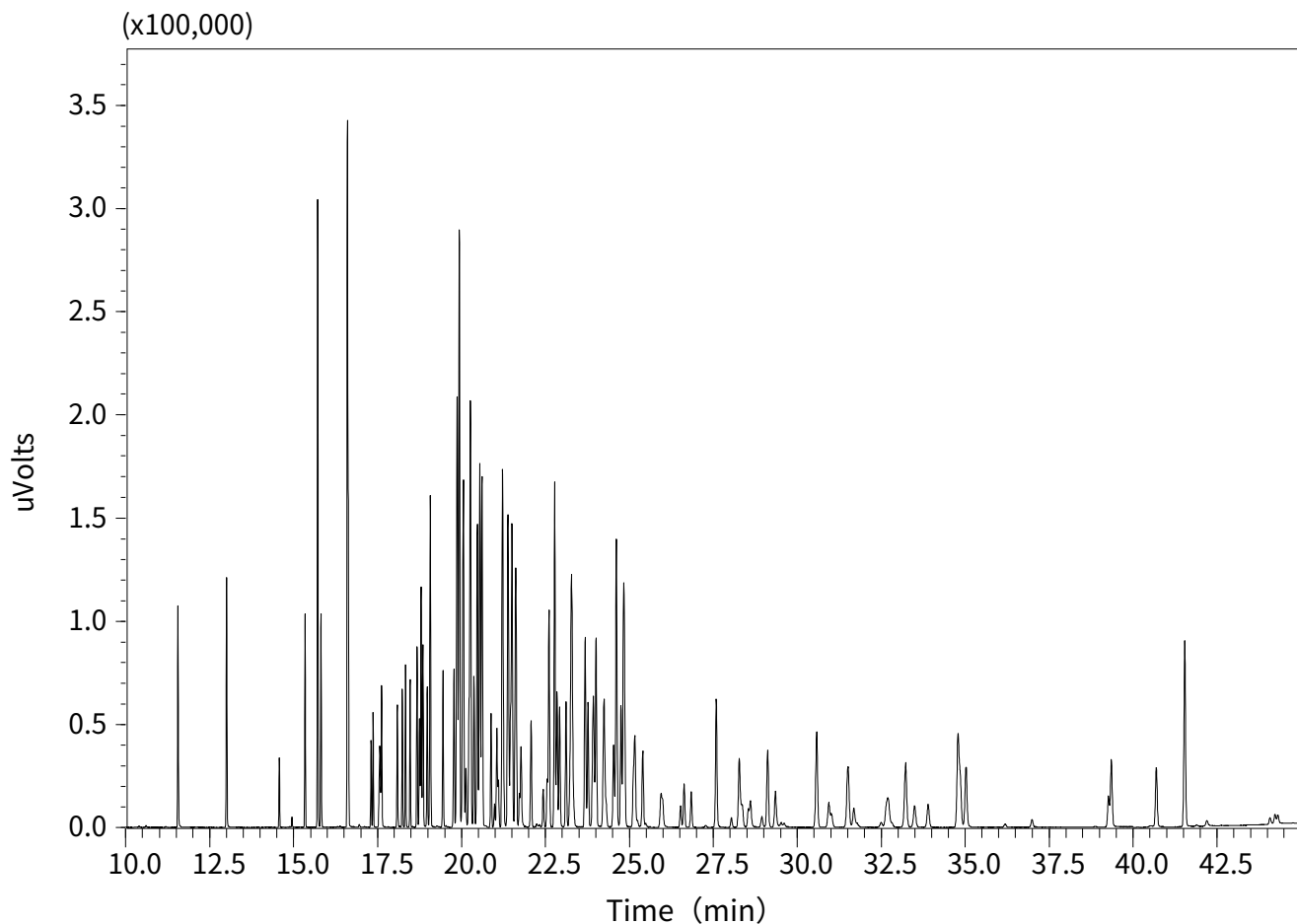

Conditions

System : GC-MS

Column : InertCap 5MS/Sil
0.25 mm I.D. x 30 m df = 0.25 μ m

Col. Cat. No. : 1010-15142

Col. Temp. : 50 °C(3 min hold) - 10 °C/min - 200 °C - 3 °C/min - 230 °C(5 min hold) -
5 °C/min - 280 °C(2 min hold)

Carrier Gas : He 100.0kPa

Injection : Splitless 1min
250 °C

Detection : MS SIM
Interface Temp. 280 °C

Sample Size : Analytes in Dichloromethane
1 μ L

InertCap Applications

No.	Compound name	RT (min)
-----	---------------	----------

Pesticides (1 $\mu\text{g/mL}$ each)