


InertSearch™ for LC

Inertsil® Applications

Analysis of Losartan Potassium (Inertsil ODS-3)

Data No. LA937-811

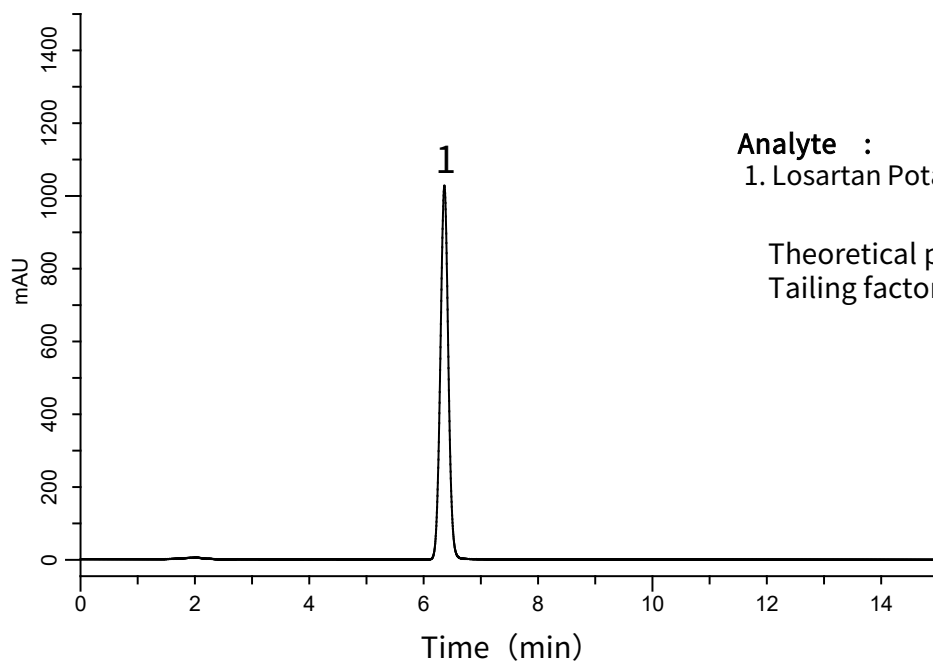

Analyte :

1. Losartan Potassium (250 mg/L)

Theoretical plates : 10,128

Tailing factor : 1.0

Conditions

System	: GL-7400 HPLC system
Column	: Inertsil ODS-3 (5 μ m, 250 x 4.0 mm I.D)
Column Cat. No.	: 5020-01762
Eluent	: A) CH ₃ CN B) 0.1% H ₃ PO ₄ A/B = 60/40, v/v
Flow Rate	: 1.0 mL/min
Col. Temp.	: 35 °C
Detection	: UV 220 nm (GL-7450 UV Detector)
Injection Vol.	: 10 μ L
Sampe	: Standard