


# InertSearch™ for LC

Inertsil® Applications

## Analysis of Spiramycins

Data No. LA991-0894

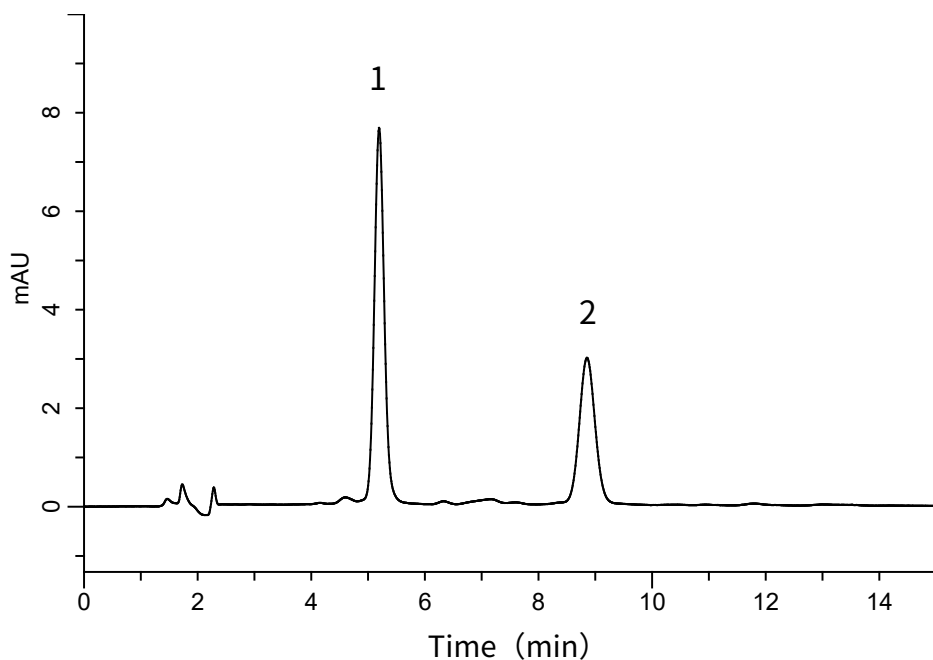

### Conditions

Column : InertSustain C18 (5 $\mu$ m, 150 x 4.6 mm I.D.)

Column Cat. No. : 5020-07345

Eluent : A) CH<sub>3</sub>CN  
B) 65 mM NaH<sub>2</sub>PO<sub>4</sub> (pH 2.5, H<sub>3</sub>PO<sub>4</sub>)  
A/B = 25/75, v/v

Flow rate : 0.5 mL/min

Col. Temp. : 40 °C

Detection : UV 235 nm

Injection Vol. : 10  $\mu$ L

Sample : Standard

### Analyte:

1. Neospiramycin I 2.0 mg/L

2. Spiramycin I 2.0 mg/L