


InertSearch™ for LC

Inertsil® Applications

Analysis of Pullulan

Data No. LB255-0871

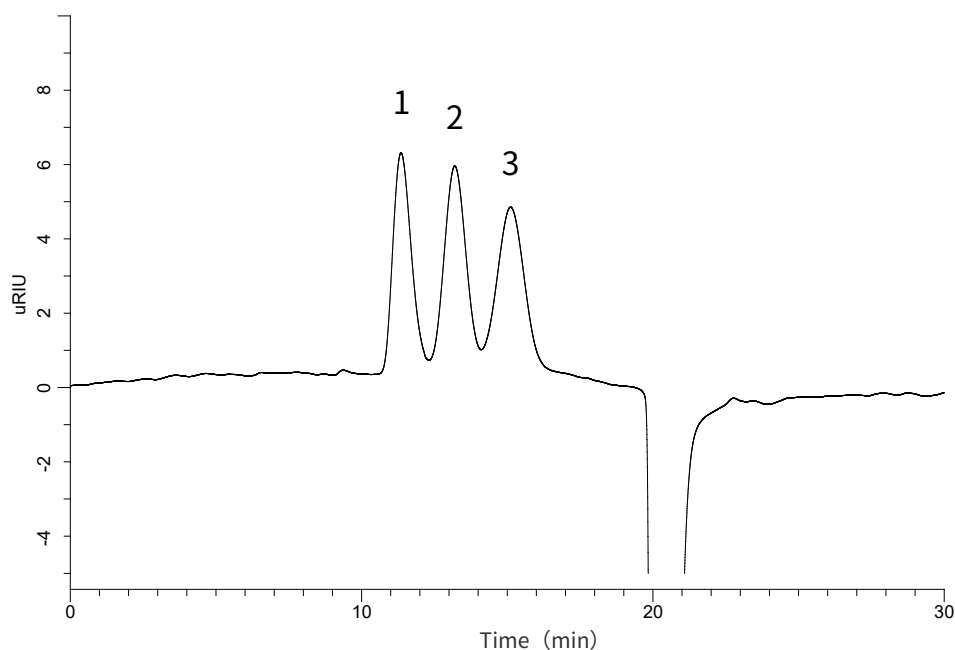

Conditions

System : GL-7400 HPLC system

Column : 1. Inertsil WP300 Diol (5 µm, 250 x 4.6 mm I.D.)
+2. Inertsil Diol (5 µm, 250 x 4.6 mm I.D.)

Column Cat. No. : 1. 5020-05946
2. 5020-05646

Eluent : A) CH₃OH
B) H₂O
A/B = 20/80, v/v

Flow Rate : 0.3 mL/min

Col. Temp. : 40 °C

Detection : RI

Injection Vol. : 20 µL

Sample : Standard

Analyte:

Pullulans :
M.W.

1. 107,000	500 µg/mL
2. 21,100	500 µg/mL
3. 5,900	500 µg/mL