


InertSearch™ for LC

Inertsil® Applications

Analysis of Ergosterol(InertSustain C18)

Data No. LB014-0811

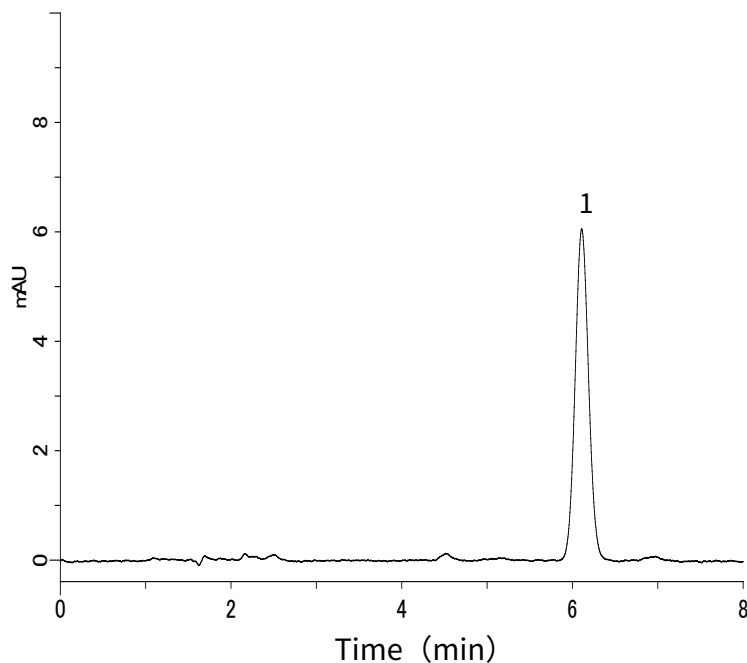

Conditions

System : GL-7400 HPLC system
Column : InertSustain C18 (5 µm, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-07345
Eluent : CH₃OH
Flow Rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 290 nm (GL-7450 UV Detector)
Injection Vol. : 10 µL
Sample : Standard

Analyte

1. Ergosterol (10 mg/L)