


Analysis of Catecholamines

Data No. LA791-0000

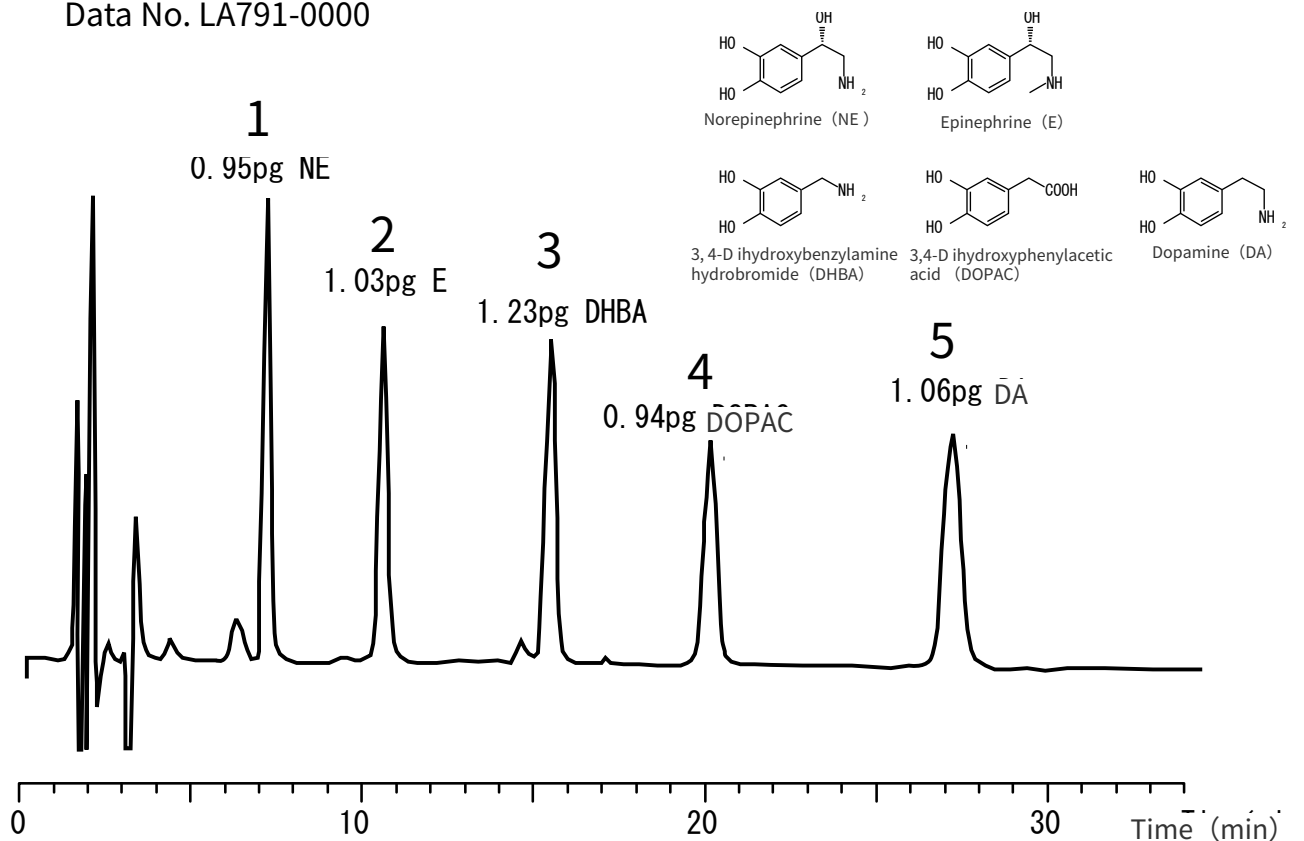

Conditions

System : AccuStream 700 series
Column : Inertsil ODS-3 (3μm, 150 x 0.3 mm I.D.)
Column Cat. No. : 5020-01871
Eluent : A) CH₃OH
 B) 50 mM Citric acid,
 50 mg/L EDTA,
 50 mg/L Octanesulfonate Na,
 50 mM H₃PO₄ (pH 3.0, KOH)
 A/B = 7/93, v/v
Flow Rate : 4.4 μL/min
Col. Temp. : Ambient
Detection : ECD 550 mV vs Ag/AgCl (Glassy Carbon)
Injection Vol. : 0.1 μL

Analyte

1. Norepinephrine (NE)
2. Epinephrine (E)
3. 3,4-Dihydroxybenzylamine hydrobromide (DHBA)
4. 3,4-Dihydroxyphenylacetic acid (DOPAC)
5. Dopamine (DA)