


InertSearch™ for LC

Inertsil® Applications

Analysis of 3 Pesticides

Data No. LA778-0000

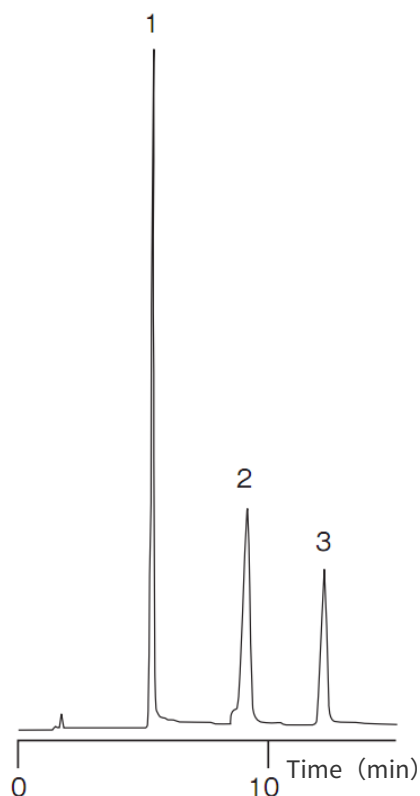

Conditions

Column : Inertsil ODS-80A (5 μ m, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-01601
Eluent : 2.0 g IPCC-MS7 in 300 mL CH₃CN + 700mL H₂O
: IPCC-MS7 : ion pair reagent for LC/MS (Cat. No. 5040-31705, 1 g)
Flow Rate : 1.0 mL/min
Col. Temp. : Ambient
Detection : UV 240 nm
Injection Vol. : 1 μ L

Analyte

1. Oxine-Copper
2. Thiabendazole (TBZ)
3. Benomyl