


InertSearch™ for LC

Inertsil® Applications

Analysis of 4 drugs for animal (animal product III)

Data No. LA827-0811

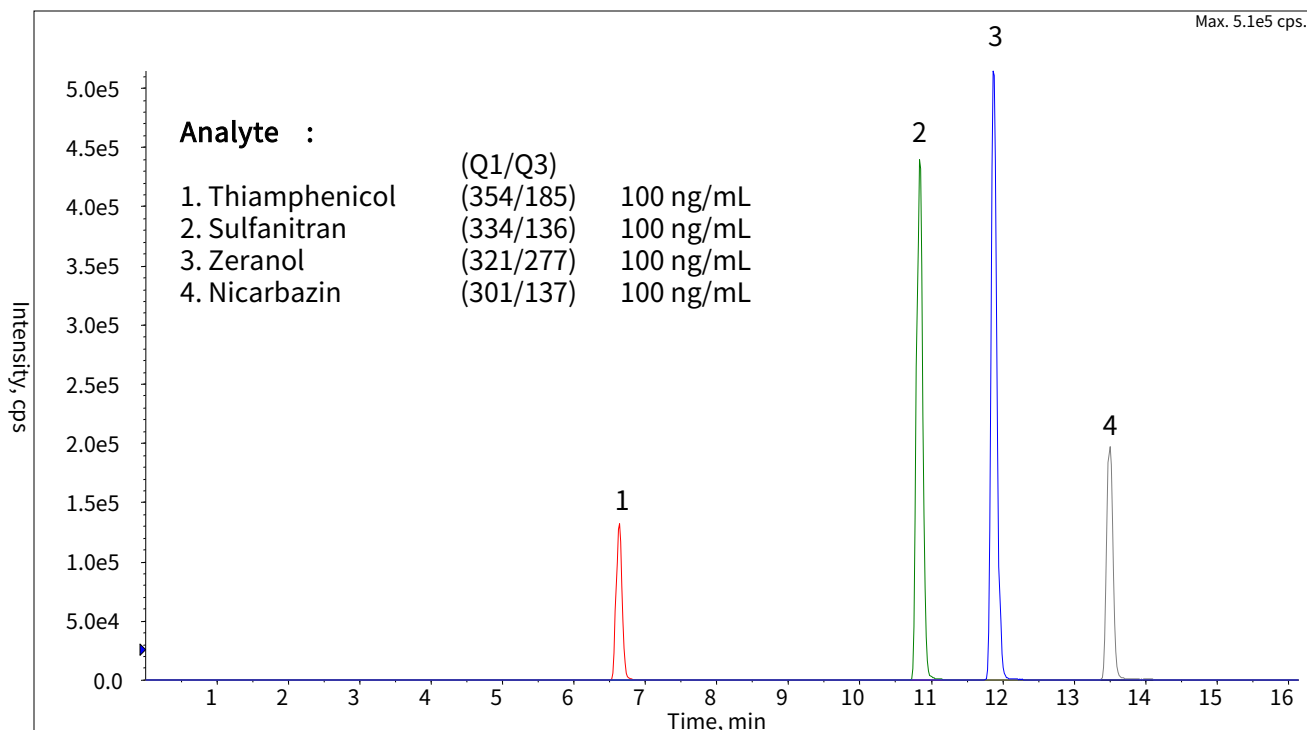

Conditions

System : LC800 HPLC system
4000 Q TRAP®

Column : Inertsil ODS-4 (3 µm, 100 x 2.1 mm I.D)

Column Cat. No. : 5020-04014

Eluent : A) 0.0025 % Formic acid in CH₃CN
B) 0.0025 % Formic acid in H₂O
A/B = 5/95 - 15 min - 80/20 - 5 min - 80/20 - 0.1 min - 5/95, v/v

Flow Rate : 0.2 mL/min

Col. Temp. : 40 °C

Detection : LC/MS/MS (4000Q TRAP®: ESI, Negative, MRM)

CUR	CAD	IS	TEM	GS1	GS2
10	5	-4500	600	70	80

Injection Vol. : 5 µL

Sample : drugs for animal

Related No. : LA820