


Analysis of Phytosterol (InertSustain C18)

Data No. LA891-0142

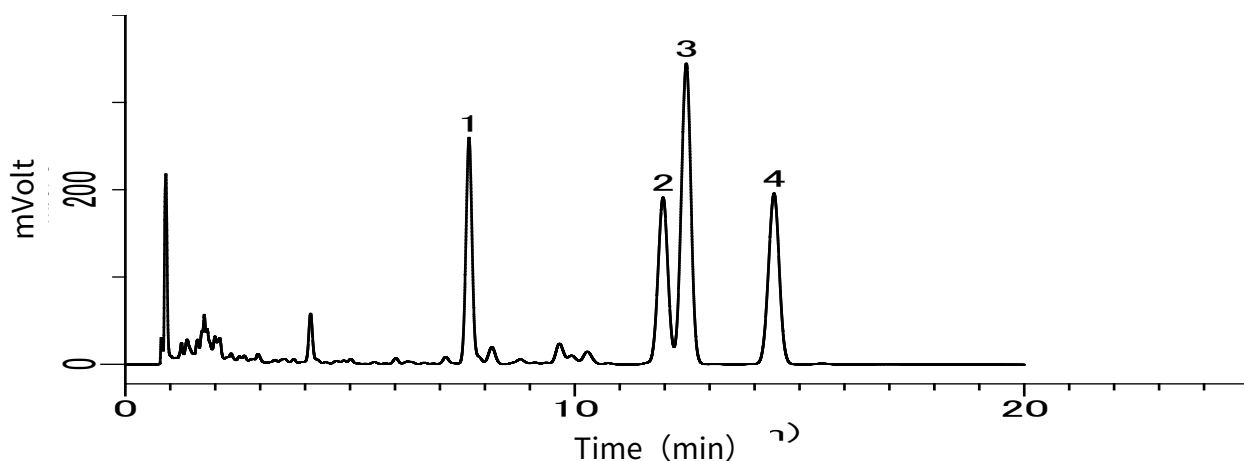

Conditions

System : LC800 HPLC system
Column : InertSustain C18 HP (3μm, 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-14415
Eluent : A) CH₃CN
 : B) H₂O
 : A/B = 95/5, v/v
Flow rate : 0.4 mL/min
Col. Temp. : 40 °C
Detection : UV 205 nm (GL-7452A PDA Detector)
Injection Vol. : 2 μL
Sample : Phytosterol


Analyte:

1. Ergosterol (1 mg/mL)
2. Campesterol (1 mg/mL)
3. Stigmasterol (1 mg/mL)
4. β-Sitosterol (1 mg/mL)


Analysis of Phytosterol (InertSustain C18)

Data No. LA891-0142


1. Ergosterol

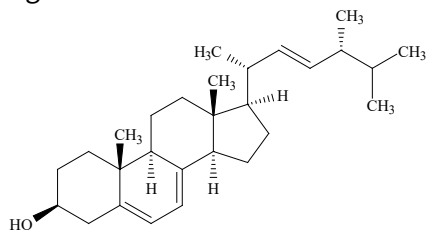

2. Campesterol

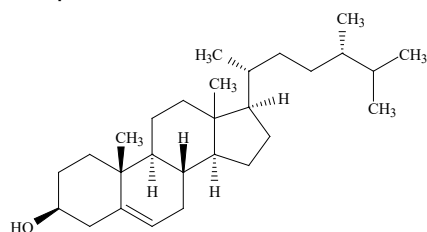

3. Stigmasterol

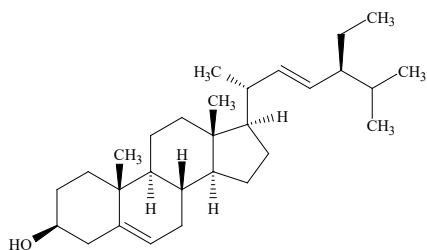

4. β -Sitosterol

