


Analysis of aliphatic acid (prelabel derivatization by ADAM reagent)

Data No. LA587-0696

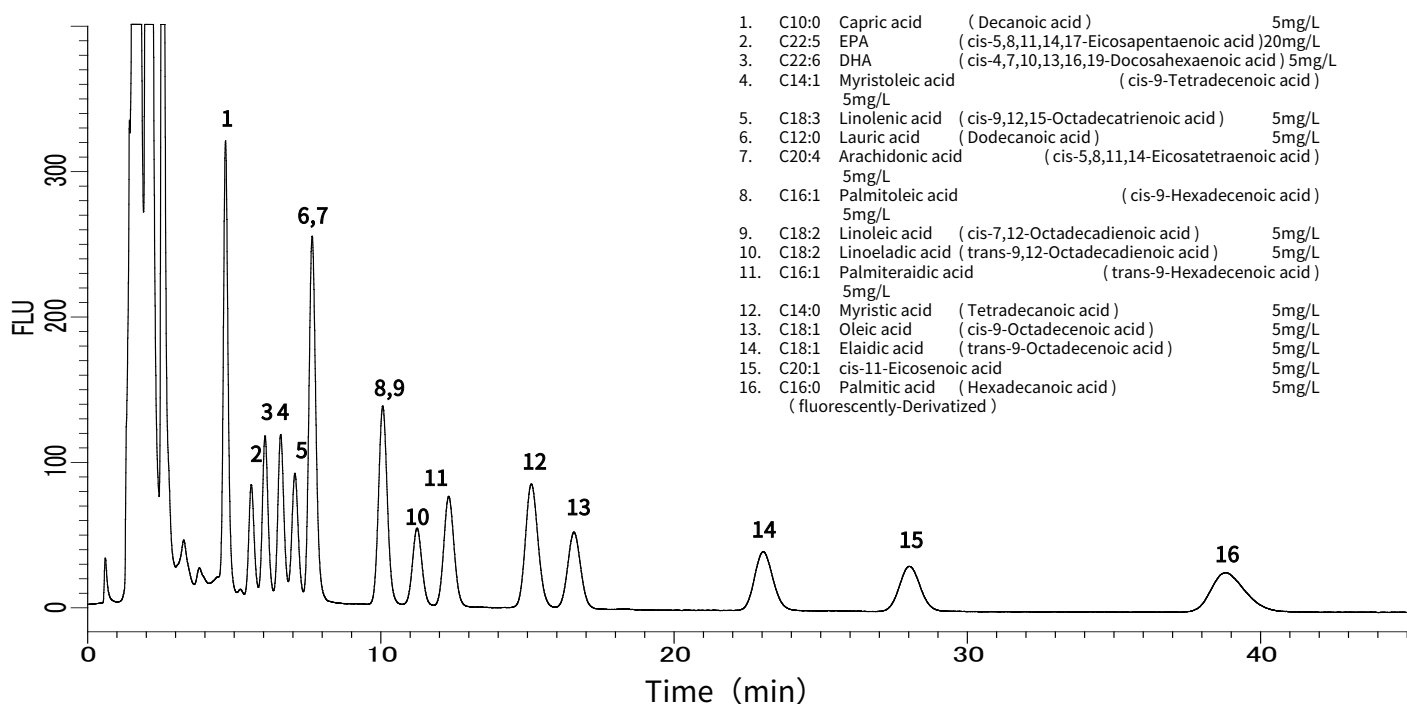

Conditions

System : GL-7400 HPLC system
Column : Inertsil ODS-P (5 µm, 250 x 3.0 mm I.D.)
Column Cat. No. :
Eluent : CH₃CN
Flow Rate : 0.9 mL/min
Col. Temp. : 20 °C
Detection : FL Ex 365 nm Em 412 nm (GL-7453A FL detector)
Injection Vol. : 5 µL
Analyte : aliphatic acid derivatized with ADAM reagent