


InertSearch™ for LC

Inertsil® Applications

Analysis of Sugar

Data No. LB140-0642

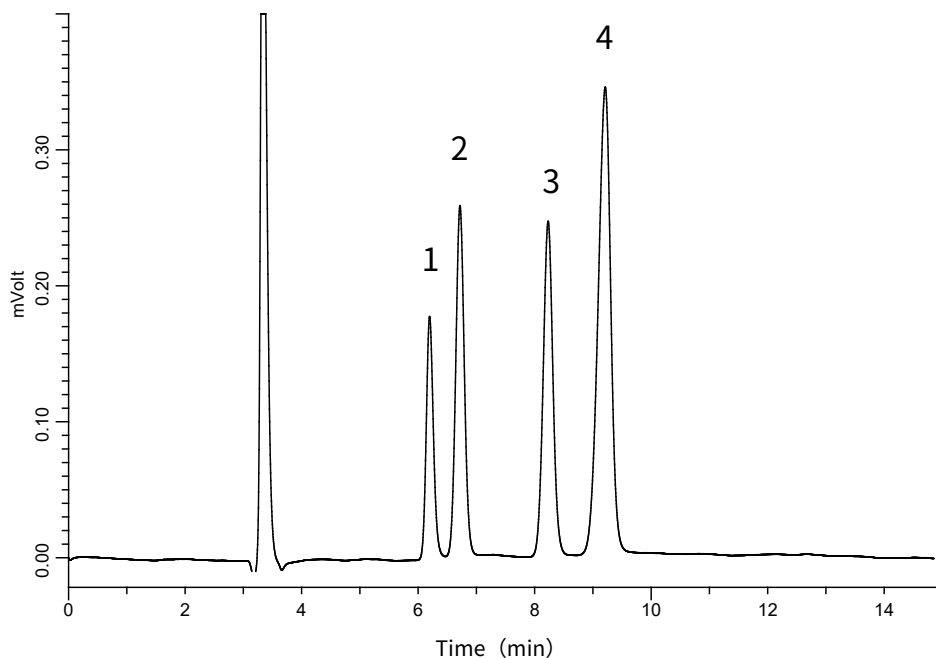

Conditions

Column : InertSustain NH₂
(5 μm, 250 x 4.6 mm I.D.)

Column Cat. No. : 5020-16628

Eluent : A) CH₃CN
B) H₂O
A/B = 75/25, v/v

Flow Rate : 1.0 mL/min

Col. Temp. : 40 °C

Detection : RI

Injection Vol. : 5 μL

Sample : Standard

Analyte:

1. Fructose
2. Glucose
3. Sucrose
4. Maltose