


Analysis of Lactosucrose

Data No. LB005-0848

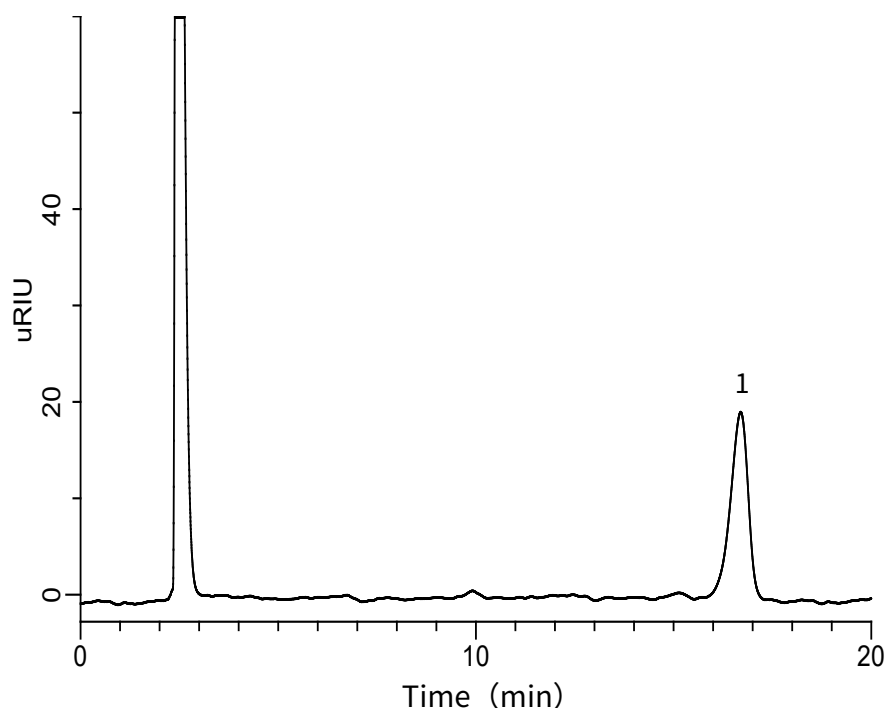

Conditions

System : GL-7400 HPLC system
Column : Inertsil Amide (5 μ m, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-07836
Eluent : A) CH₃CN
 : B) H₂O
 : A/B = 71/29, v/v
Flow rate : 1.2 mL/min
Col. Temp. : 35 °C
Detection : RI (GL-7454 RI Detector)
Injection Vol. : 20 μ L
Sample : Standard

Analyte:

1. Lactosucrose 5 mg/mL

