


InertSearch™ for LC

Inertsil® Applications

Analysis of Pitavastatin Calcium

Data No. LB027-0894

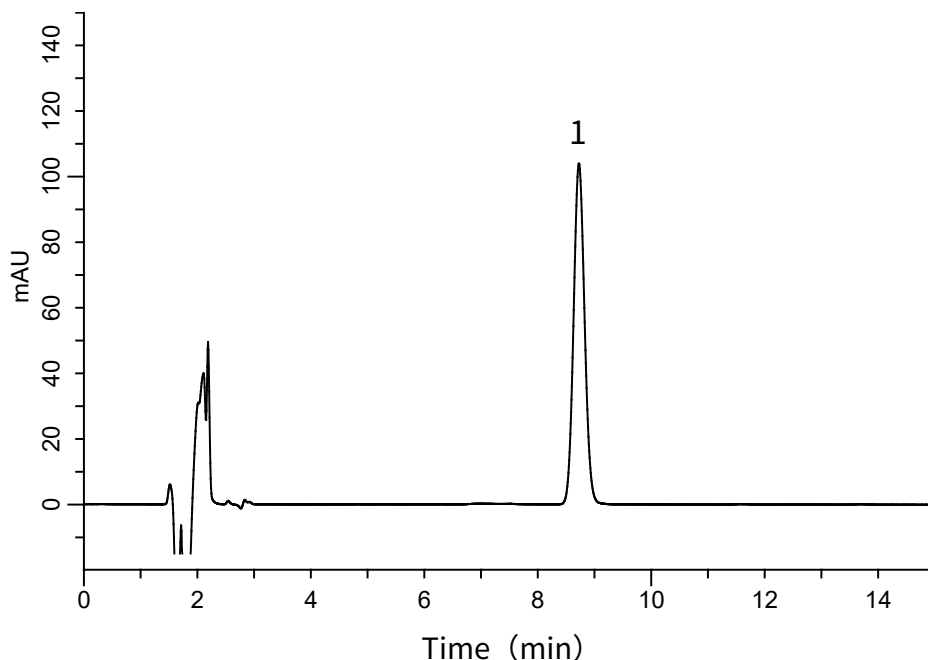

Conditions

System : GL-7400 HPLC system
Column : InertSustain C18 (5μm, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-07346
Eluent : A) THF
 : B) CH₃CN
 : C) Citric acid buffer*
 : A/B/C = 20/27/53, v/v/v (Premix)
Flow rate : 1.2 mL/min
Col. Temp. : 30 °C
Detection : UV 244 nm (GL-7452A PDA Detector)
Injection Vol. : 50 μL

Analyte:

1. Pitavastatin calcium 0.006 mg/mL

* Dissolve 10.5 g of citric acid monohydrate in 900 mL of water.
Adjust pH 4.0 by 28 % ammonia solution.
Add water to make 1000 mL.