


InertSearch™ for LC

Inertsil® Applications

Analysis of Estrogen

Data No. LB137-0642

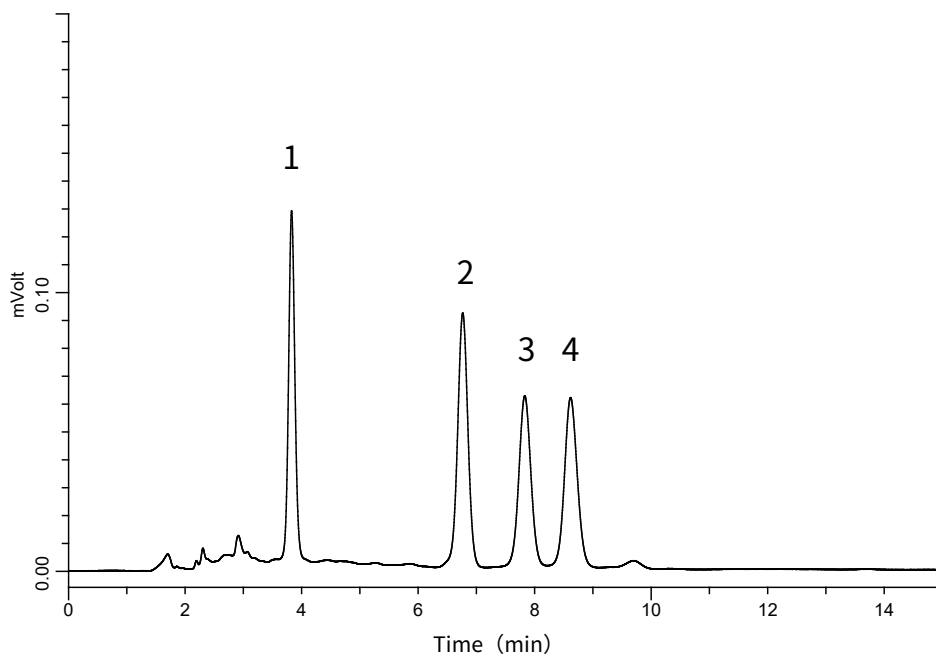

Conditions

Column : InertSustain Phenyl
(5 µm, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-16327
Eluent : A) CH₃OH
B) H₂O
A/B = 60/40, v/v
Flow Rate : 0.8 mL/min
Col. Temp. : 40 °C
Detection : UV 220 nm
Injection Vol. : 5 µL
Sample : Standard

Analyte:

1. Estriol
2. Estradiol
3. Ethynylestradiol
4. Estrone