


InertSearch™ for LC

Inertsil® Applications

Analysis of Sugar Alcohols

Data No. LA926-0812

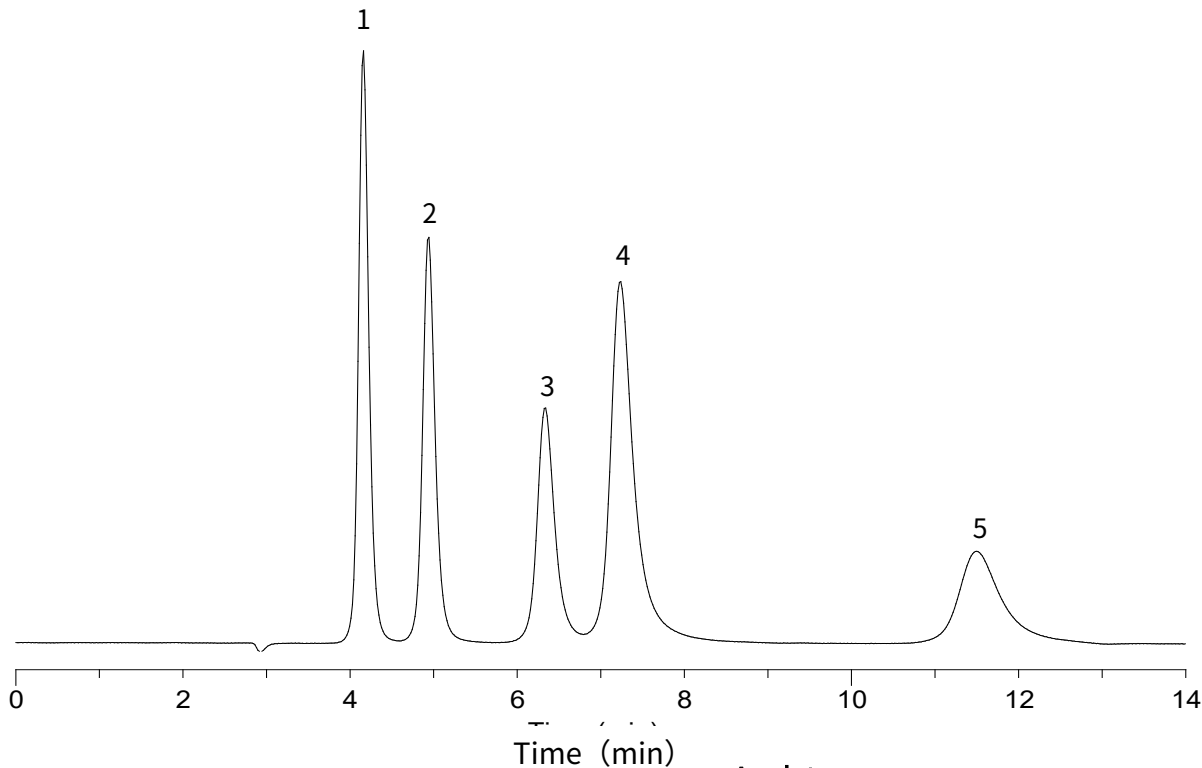

Conditions

System : LC700 system
Column : InertSphere Sugar-1 (5µm, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-11001
Eluent : 100 mM NaOH
Flow rate : 0.5 mL/min
Col. Temp. : 25°C
Detection : ECD Pulse Mode (ED723, Gold)

E1:	150 mV	t1:	600 ms
E2:	-1500 mV	t2:	50 ms
E3:	600 mV	t3:	50 ms
E4:	-200 mV	t4:	100 ms
		ts:	50 ms

Analyte:

1. Inositol	(10 mg/L)
2. Erythritol	(10 mg/L)
3. Arabitol	(10 mg/L)
4. Galactitol	(10 mg/L)
5. Lactitol	(10 mg/L)

Injection Vol. : 10 µL
Sample : Sugar alcohols