


Analysis of Peptides (MonoCap C18 for Nano-flow)

Data No. LA797-0700

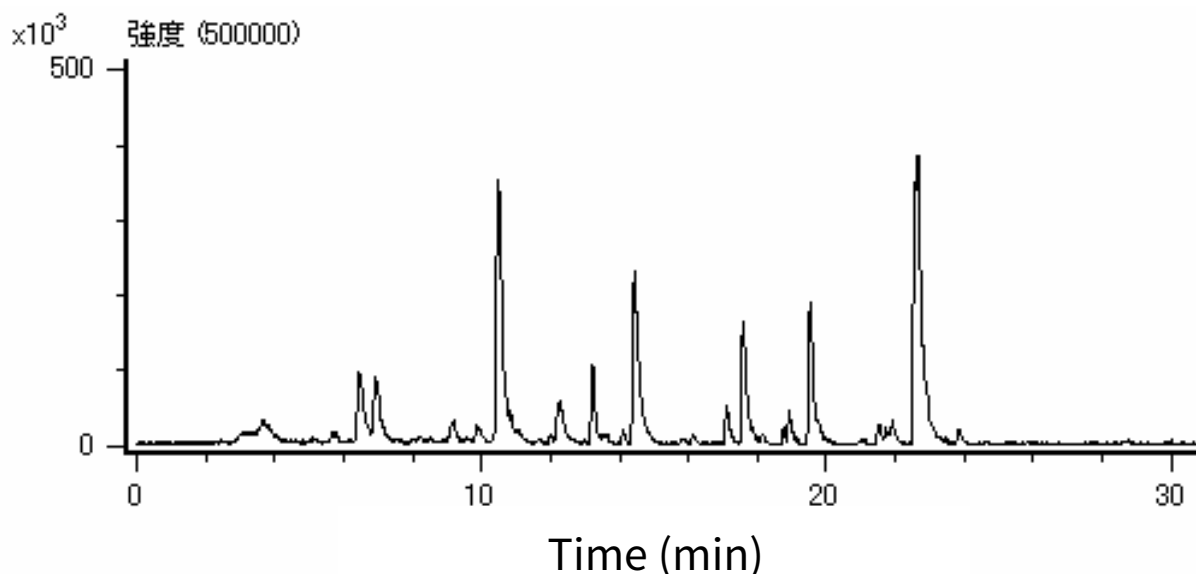

Conditions

System : AccuStream 700 series
Column : MonoCap C18 for Nano-flow (50 x 0.1 mm I.D.)
Column Cat. No. : 5020-10153
Sprayer : MonoSpray FS (50 mm x 50 μ m I.D.)
Cat. No. : 5010-20006
Eluent : A) 0.1 % HCOOH in CH₃CN
B) 0.1 % HCOOH in H₂O
A/B = 10/90 – 40 min – 50/50, v/v
Flow Rate : 0.5 μ L/min
Col. Temp. : Ambient
Detection : Nano-ESI MS positive TIC (JMS-T100-LC, m/z: 400 – 1400)
Injection Vol. : 0.1 μ L

Analyte

Lys Endopeptidase Digest of
 β -Lactoglobulin