


Analysis of Fatty acid methyl ester

Data No. LB169-0521

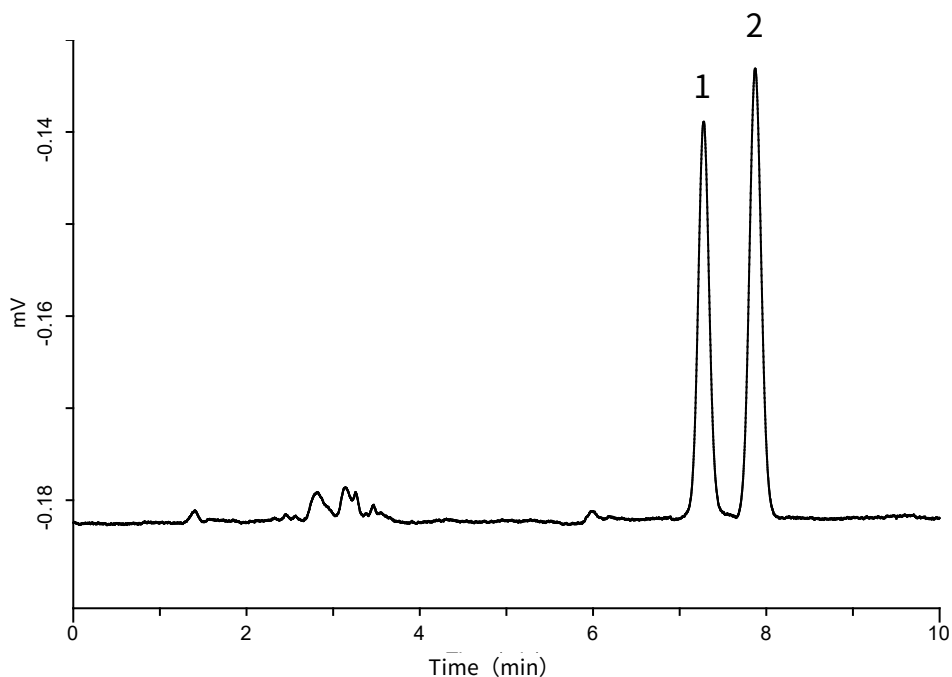

Conditions

System : GL-7400 HPLC system
Column : Inertsil C30 S-Select
(5 μ m, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-17037
Eluent : CH₃OH
Flow Rate : 1.0 mL/min
Col. Temp. : 35 °C
Detection : UV 205 nm (GL-7452 PDA Detector)
Injection Vol. : 10 μ L
Sample : Standard

Analyte:

1. Oleic acid methyl ester	500 mg/L
2. Elaidic acid methyl ester	500 mg/L