


Analysis of Calcium pantothenate and Biotin

Data No. LA884-0142

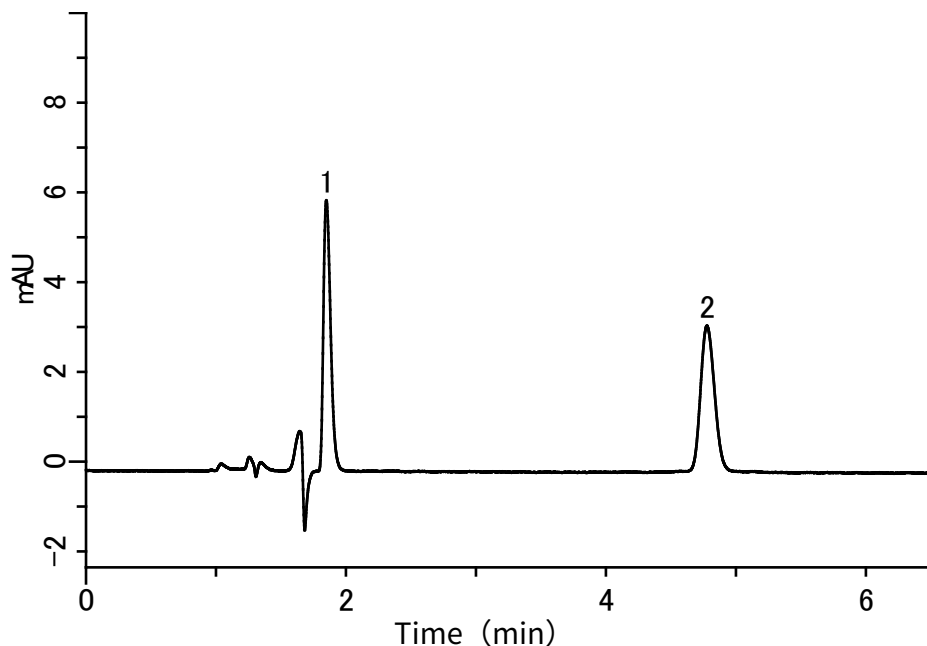

Conditions

System : GL-7400 HPLC system
Column : Inertsil ODS-4 HP (3μm, 150 x 3.0 mm I.D.)
Column Cat. No. : 5020-14005
Eluent : A) CH₃CN
 : B) 0.1% H₃PO₄
 : A/B = 15/85, v/v
Flow rate : 0.6 mL/min
Col. Temp. : 40 °C
Detection : UV 210 nm (GL-7452A PDA Detector)
Injection Vol. : 10 μL
Sample : Standard

Analyte:

1. Calcium pantothenate (10 μg/mL)
2. Biotin (10 μg/mL)