


InertSearch™ for LC

Inertsil® Applications

Analysis of Acetic acid and Acrylamide

Data No. LB266-0659

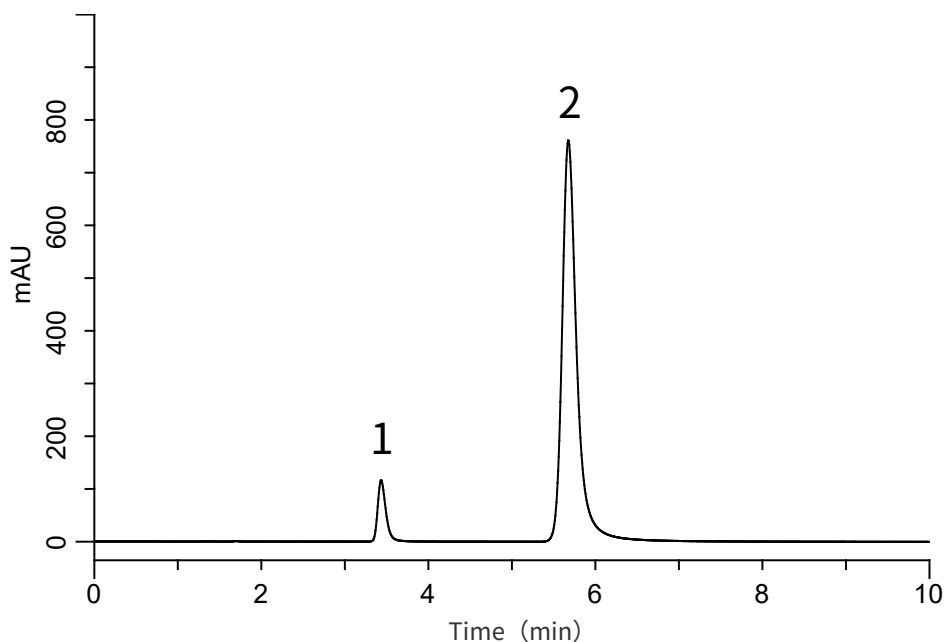

Conditions

System : GL-7400 HPLC system
Column : Inertsil C30 S-Select
(5 μ m, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-17036
Eluent : 0.1 % H_3PO_4 in H_2O
Flow Rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 210 nm (GL-7452A PDA Detector)
Injection Vol. : 5 μ L
Sample : Standard

Analyte:

1. Acetic acid	5.0 μ L/mL
2. Acrylamide	0.5 mg/mL