


InertSearch™ for LC

Inertsil® Applications

Analysis of Catecholamines

Data No. LA857-0871

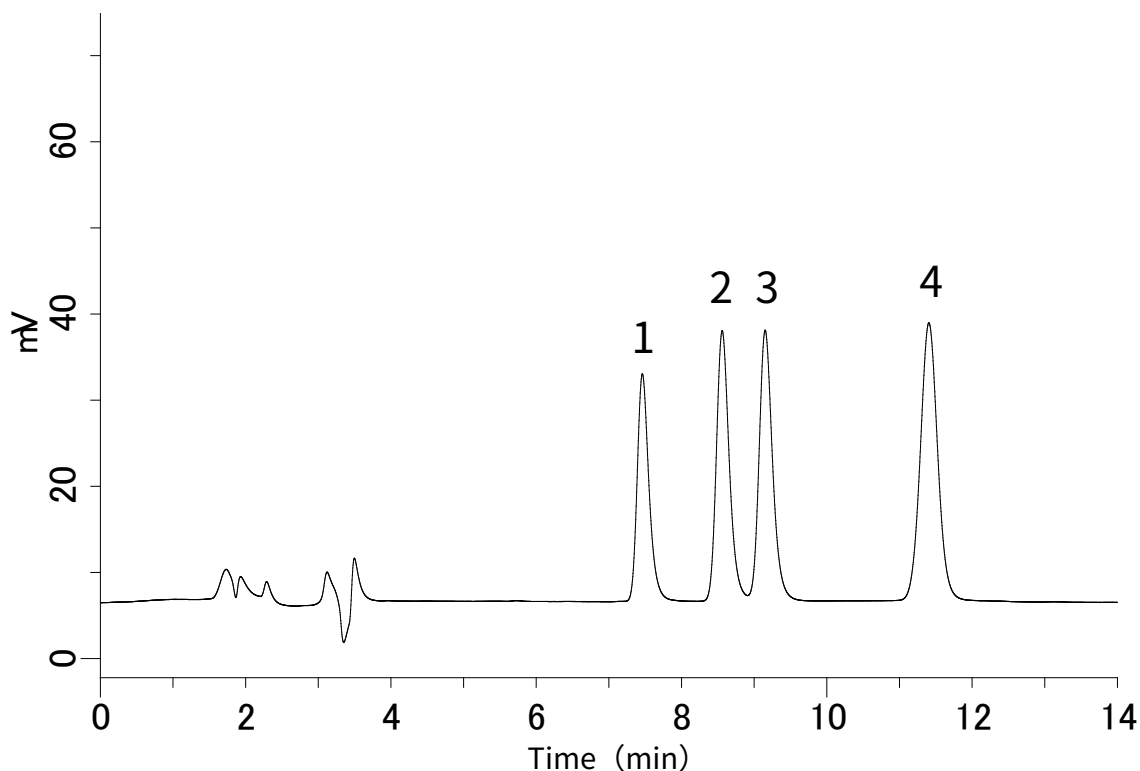

Conditions

System : LC700 HPLC system
Column : Inertsil ODS-4 (5µm, 250 x 3.0 mm I.D.)
Column Cat. No. : 5020-03926
Eluent : A) CH₃CN
 : B) 20 mM Sodium acetate + 20 mM citric acid
 : + 1 g/L IPCC-08 in water
 : (IPCC-08: Sodium 1-Octane sulfonate)
 : A/B = 16/100, v/v (Premix)

Flow rate : 0.5 mL/min

Col. Temp. : 35°C

Detection : ECD 800 mV vs. Ag/AgCl (ED703 ECDetector, Diamond)

Injection Vol. : 20 µL

Sample : Catecholamines

Analyte:

1. Norepinephrine	(100 µg/L)
2. Epinephrine	(100 µg/L)
3. 3,4-Dihydroxybenzylamine	(100 µg/L)
4. Dopamine	(100 µg/L)