


InertSearch™ for LC

Inertsil® Applications

Analysis of Roxithromycin

(Under the Condition of the Japanese Pharmacopoeia)

Data No. LB329-0919

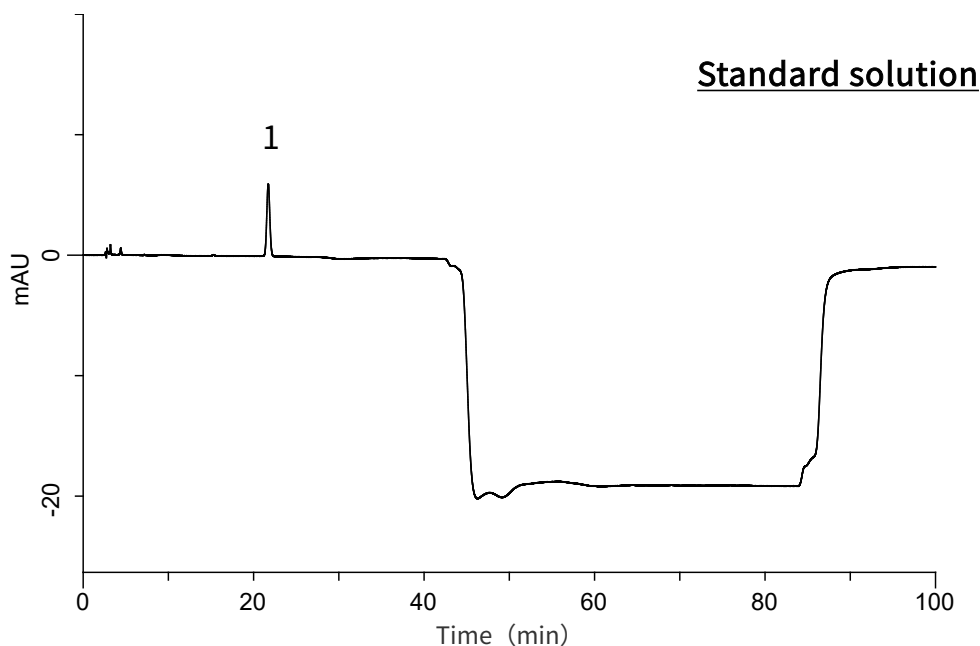

Conditions

System : GL7700 HPLC system
Column : InertSustainSwift C18
(5 μ m, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-88027
Eluent : A) Buffer*
B) CH₃CN/H₂O = 70/30, v/v

Time (min)	A (vol%)	B (vol%)
0.0	100	0
38.0	100	0
39.0	90	10
80.0	90	10
80.1	100	0
100.0	100	0

Flow Rate : 0.92 mL/min
Col. Temp. : 25 °C
Detection : UV 205 nm (UV7750 UV Detector)
Injection Vol. : 20 μ L
Sample : Standard

*Dissolve 34 g of ammonium dihydrogenphosphate in 710 mL of water.
Adjust pH 5.3 by 2 mol/L sodium hydroxide test solution.
Add 315 mL of acetonitrile.

Analyte:

1. Roxithromycin 20 mg/L

Theoretical plates : 17,431 ($\geq 9,000$)
Tailing factor : 1.05 (≤ 1.5)
RSD of the
peak area (%) (n=5) : 0.09 (≤ 2.0)

【NOTE】

- 1) The retention time will be shifted easily under this conditions.
- 2) Salting-out will be occurred due to high concentration buffer.