


InertSearch™ for LC

Inertsil® Applications

Analysis of Vitamin D

Data No. LB467-0812

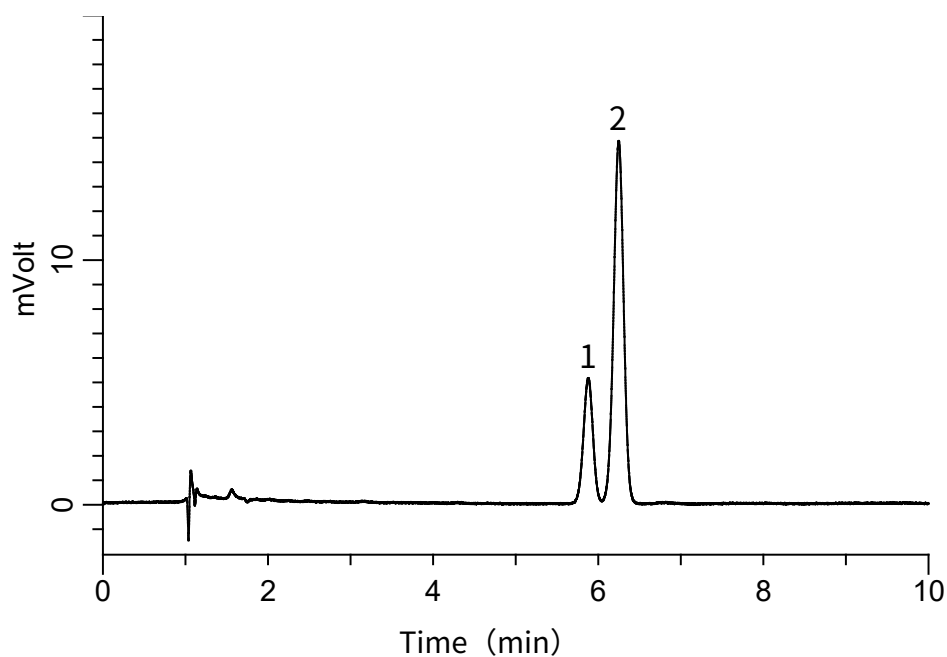

Conditions

Column : Inertsil ODS-HL
(3 μ m, 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-87231
Eluent : CH₃OH
Flow rate : 0.3 mL/min
Col. Temp. : 25 °C
Detection : UV 265 nm
Injection Vol. : 5 μ L
Sample : Standard

Analyte:

1. Vitamin D2
2. Vitamin D3
(5 mg/L each)