


Analysis of Chelating compounds

Data No. LB340-0806

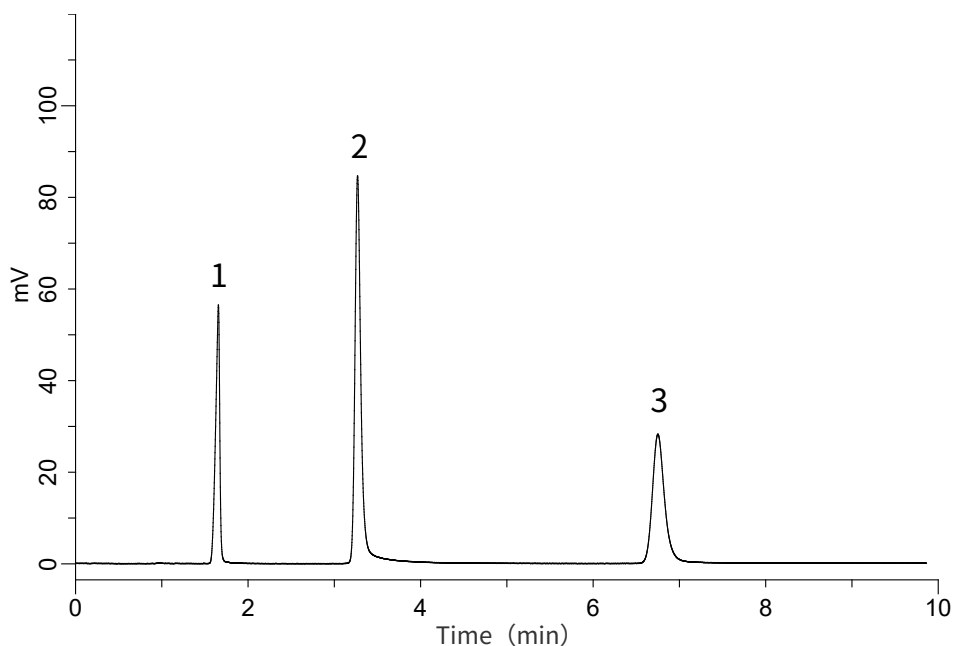


Conditions


Column : InertSustain C18
(5 µm, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-07345
Eluent : A) CH₃CN
B) 0.1 % HCOOH in H₂O
A/B = 60/40, v/v
Flow Rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 310 nm
Injection Vol. : 2 µL
Sample : Standard

Analyte:


1. 8-Quinolinol
2. Hinokitiol
3. Piroctone Olamine
(100 µg/mL each)

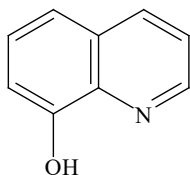

1. 8-Quinolinol

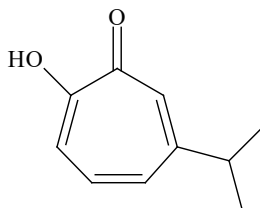

2. Hinokitiol

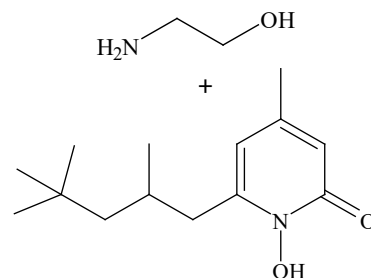

3. Piroctone Olamine