


InertSearch™ for LC

Inertsil® Applications

Analysis of Acesulfame potassium

(Under the Condition of the JPE (Japanese Pharmaceutical Excipients))

Data No. LB312-0894

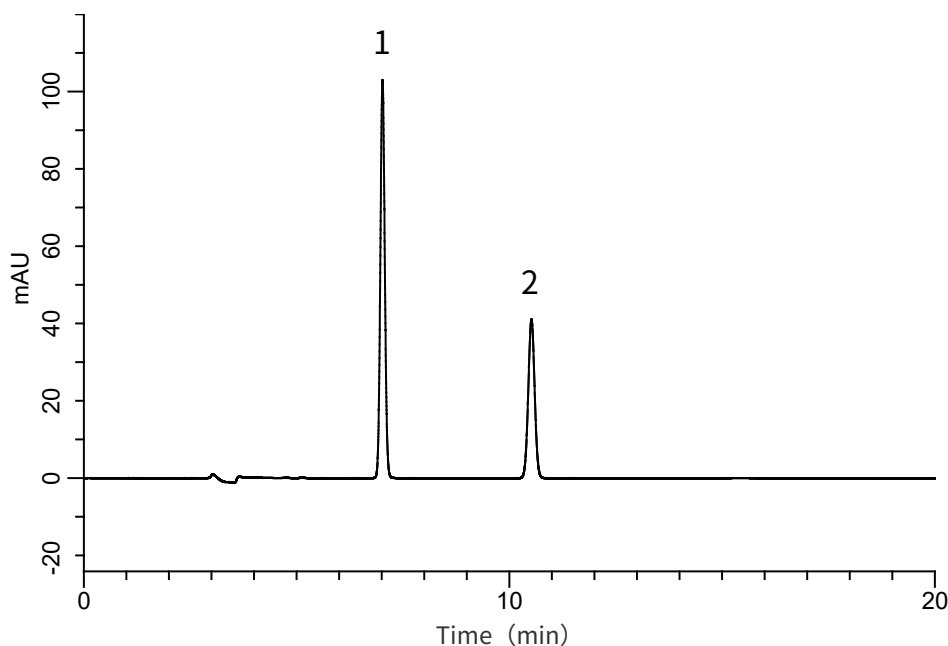

Conditions

System : GL7700 HPLC system

Column : InertSustain C18
5 μ m, 250 x 4.6 mm I.D.)

Column Cat. No. : 5020-07346

Eluent : A) CH₃CN
B) 10 mM Tetrabutylammonium hydrogensulfate
A/B = 2/3, v/v

Flow Rate : 0.83 mL/min

Col. Temp. : 40 °C

Detection : UV 227 nm (UV7750 UV Detector)

Injection Vol. : 20 μ L

Sample : Standard

Analyte:

1. Acesulfame potassium

10 mg/L

2. *p*-Hydroxybenzoic acid ethyl ester

20 mg/L

Resolution (1, 2) : 15.0 ($\geq$ 4)