


Analysis of Nucleosides

Data No. LB490-806

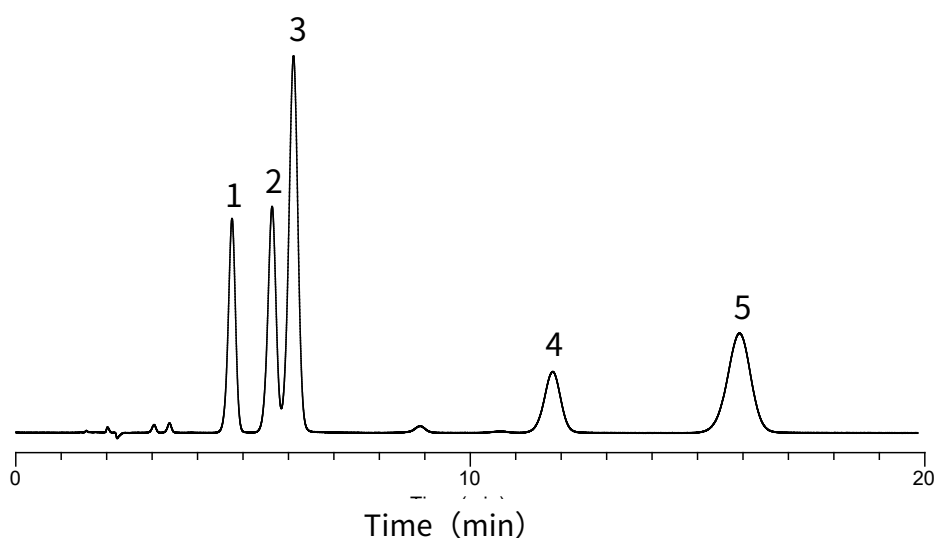

Conditions

Column : InertSustain Amide
(5 µm, 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-88607
Eluent : A) CH₃CN
B) 0.1 % HCOOH, 20 mM HCOONH₄ in H₂O
A/B = 90/10, v/v
Flow Rate : 0.2 mL/min
Col. Temp. : 40 °C
Detection : UV 254 nm
Injection Vol. : 1 µL
Sample : Standard

Analyte:

1. 5-Methyluridine
2. Uridine
3. Adenosine
4. Cytidine
5. Guanosine
(100 mg/L each)