

InertSearch™ for LC


Inertsil® Applications

Analysis of Candesartan cilexetil

(Under the Condition of the Japanese Pharmacopoeia,
Candesartan cilexetil and Hydrochlorothiazide Tablets)

Data No. LB302-0894

System performance

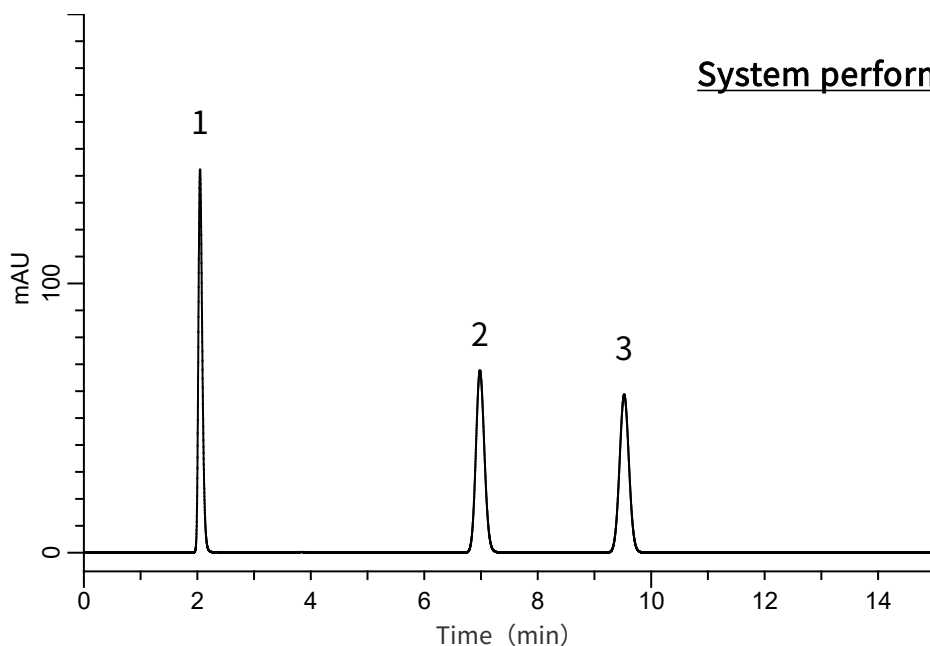

Conditions

System : GL7700 HPLC system
Column : Inertsil ODS-4
(4 μ m, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-89608
Eluent : A) CH₃CN
B) Buffer*
A/B = 11/9, v/v
Flow Rate : 0.98 mL/min
Col. Temp. : 25 °C
Detection : UV 254 nm (UV7750 UV Detector)
Injection Vol. : 10 μ L
Sample : Standard

*Dissolve 7.80 g of sodium dihydrogenphosphate dihydrate in 900 mL of water.
Adjust pH 5.5 by sodium hydroxide.
Add water to make 1,000 mL.

Analyte:

1. Hydrochlorothiazide	62 mg/L
2. Candesartan cilexetil	40 mg/L
3. Benzophenone	10 mg/L

Resolution (1, 2) : 24.3 (≥ 7)

Resolution (2, 3) : 8.6 (≥ 6)