


InertSearch™ for LC

Inertsil® Applications

Analysis of Amygdalin

Data No. LB332-0848

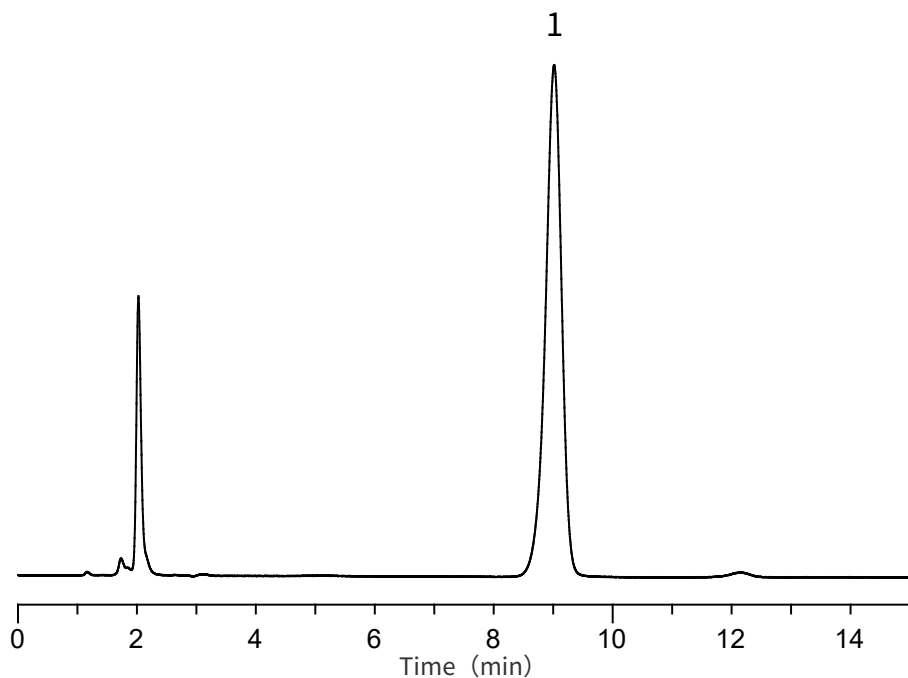

Conditions

System : GL7700 HPLC system
Column : Inertsil Amide
(5 µm, 150 x 4.6 mm I.D.)
Column Cat. No. : 5020-07835
Eluent : A) CH₃CN
B) H₂O
A/B = 88/12, v/v
Flow Rate : 1.0 mL/min
Col. Temp. : 40 °C
Detection : UV 210 nm (UV7750 UV Detector)
Injection Vol. : 10 µL
Sample : Standard

Analyte:

1. Amygdalin 50 mg/L