


Analysis of Tasty Components

Data No. LB498-806

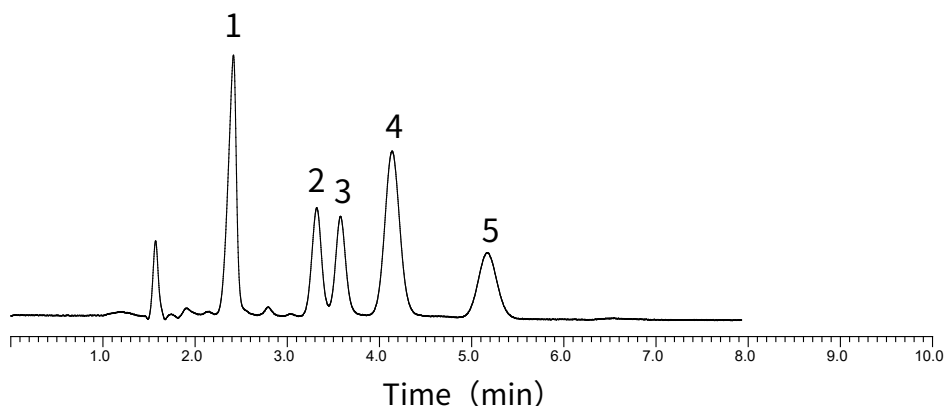

Conditions

Column : InertSustain Amide
(5 μ m, 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-88607
Eluent : A) CH₃CN
B) 25 mM Phosphate Buffer
(K₂HPO₄ and KH₂PO₄)
A/B = 75/25, v/v
Flow Rate : 0.2 mL/min
Col. Temp. : 40 °C
Detection : UV 210 nm
Injection Vol. : 1 μ L
Sample : Standard

Analyte:

1. Succinic acid	5 mg/mL
2. Glutamic acid	5 mg/mL
3. Aspartic acid	5 mg/mL
4. Inosinic acid (IMP)	0.1 mg/mL
5. Guanylic acid (GMP)	0.1 mg/mL