


Analysis of Chlorophenols

Data No. LB609-0806

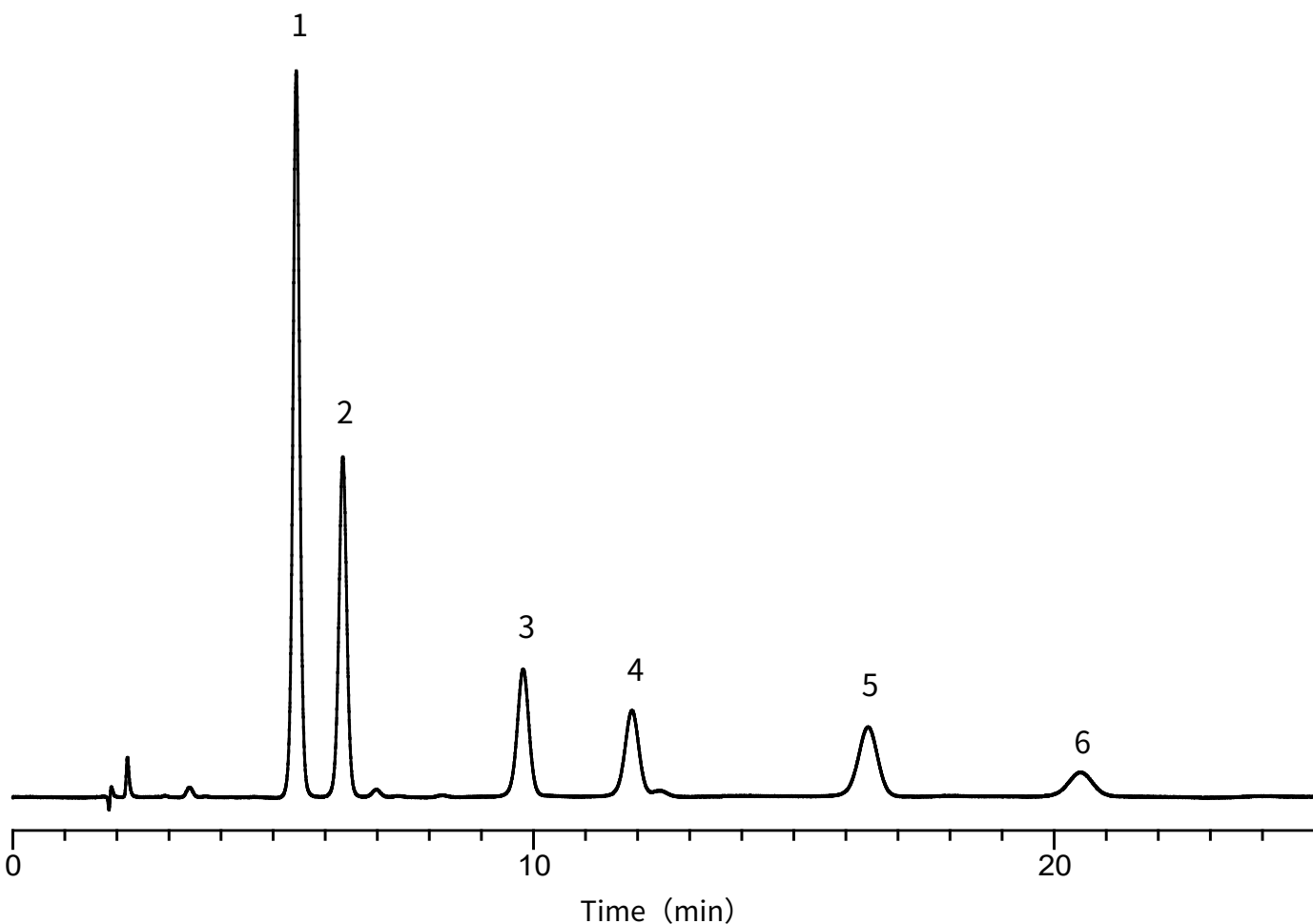

Conditions

Column : InertSustain Cyano (5 µm, 150 x 4.6 mm I.D.)

Column Cat. No. : 5020-89276

Eluent : A) CH₃OH
B) 0.1% H₃PO₄ in H₂O
A/B = 40/60, v/v

Flow rate : 1.0 mL/min

Col. Temp. : 40 °C

Detection : UV 280 nm

Injection Vol. : 2 µL

Sample : Standard

Analyte:

1. 2,6-Dichlorophenol 100 mg/L

2. 2,4-Dichlorophenol 100 mg/L

3. 2,4,5-Trichlorophenol 100 mg/L

4. 3,4,5-Trichlorophenol 100 mg/L

5. 2,3,4,5-Tetrachlorophenol 100 mg/L

6. Pentachlorophenol 100 mg/L