


Analysis of Corticosteroids

Data No. LB494-812

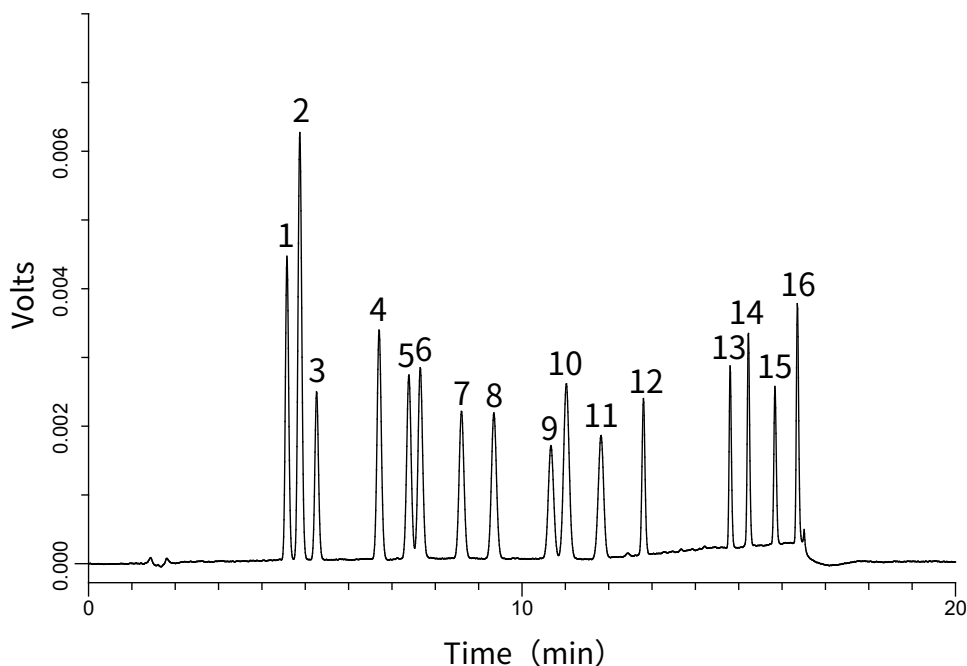

Conditions

Column : Inertsil ODS-HL
(3 µm, 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-87231
Eluent : A) CH₃CN
B) H₂O

Time (min)	A (%)	B (%)
0.00	35	65
5.00	38	62
10.0	42	58
15.0	90	10
15.1	35	65
25.1	35	65

Flow Rate : 0.3 mL/min
Col. Temp. : 30 °C
Detection : UV 260 nm
Injection Vol. : 1 µL
Sample : Standard

Analyte:

1. Prednisolone
2. Hydrocortisone
3. Cortisone
4. 6- α -Methylprednisolone
5. Betamethasone
6. Dexamethasone
7. Corticosterone
8. Triamcinolone Acetonide
9. Fluocinolone Acetonide
10. Prednisolone Acetate
11. Hydrocortisone Acetate
12. Cortisone Acetate
13. Fluocinonide
14. Betamethasone Valerate
15. Deoxycorticosterone
16. Beclomethasone