


InertSearch™ for LC

Inertsil® Applicationsv

Analysis of Eugenol

Data No. LB587-0806

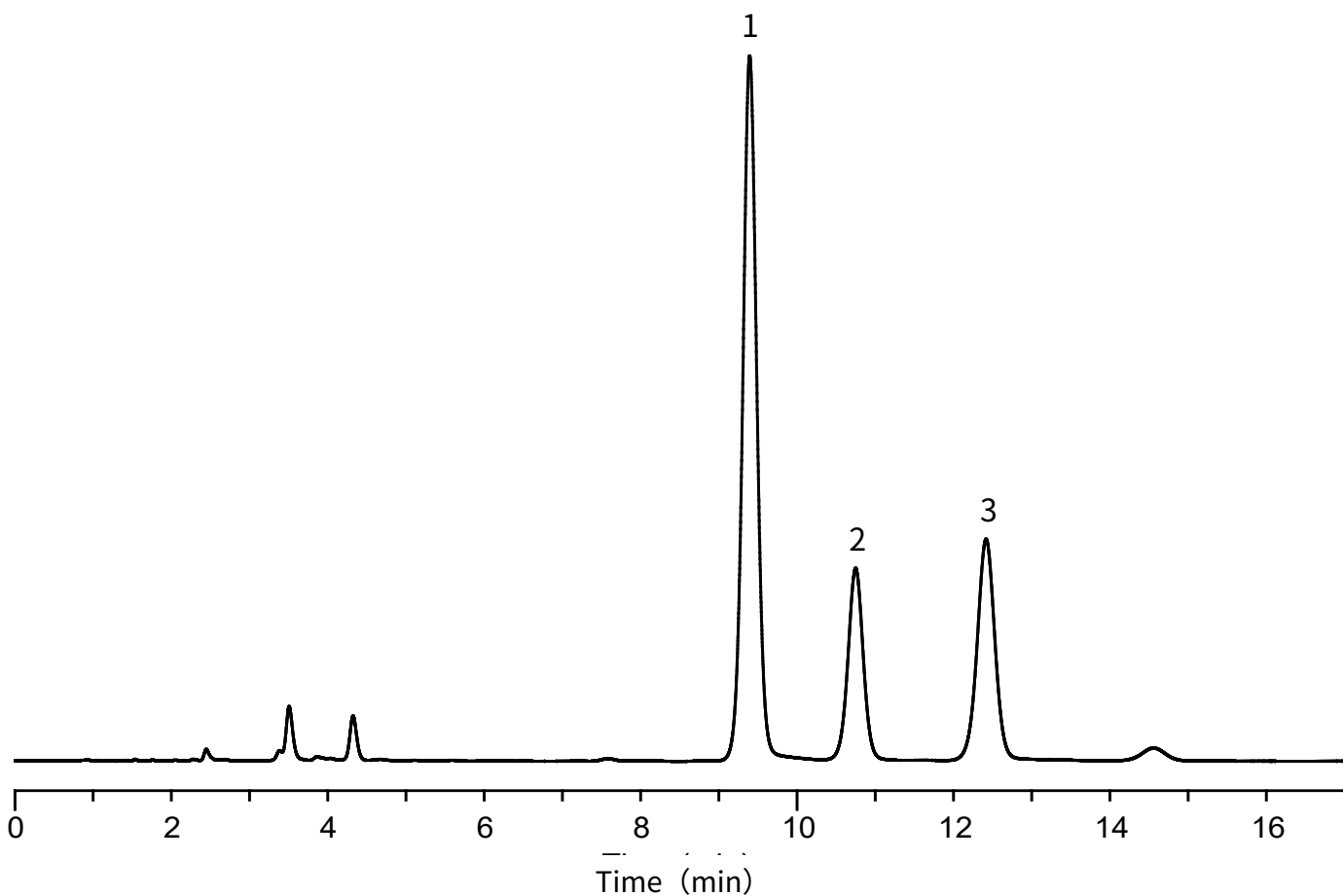

Conditions

System : GL7700 HPLC system

Column : InertSustain PFP HP (3 µm, 150 x 2.1 mm I.D.)

Column Cat. No. : 5020-87921

Eluent : A) CH₃OH

B) H₂O

A/B = 40/60, v/v

Flow rate : 0.2 mL/min

Col. Temp. : 40 °C

Detection : UV 210 nm (UV7750 with semi-micro flow cell)

Injection Vol. : 2.5 µL

Sample : Standard

Analyte:

1. Eugenol 100 mg/L

2. *cis*-Isoeugenol 100 mg/L

3. *trans*-Isoeugenol 100 mg/L