



Analysis of stable radical DPPH

Here we introduce the analysis of the 1,1-diphenyl-2-picrylhydrazyl free radical (DPPH radical), which is known as a stable radical. When analyzed using reversed-phase HPLC, a peak should be 1,1-diphenyl-2-picrylhydrazine (DPPH) was detected as a DPPH radical impurity. For this reason, we checked the retention time and the UV-VIS spectrum obtained with a diode array detector. Acetonitrile dehydrated with molecular sieves was used for sample dilution. From the calibration curve and the chromatogram at low concentration, a certain amount was decomposed before injection, and it was confirmed that it changed even in the HPLC system. (K. Suzuki)

Analysis example of DPPH free radical standard solution



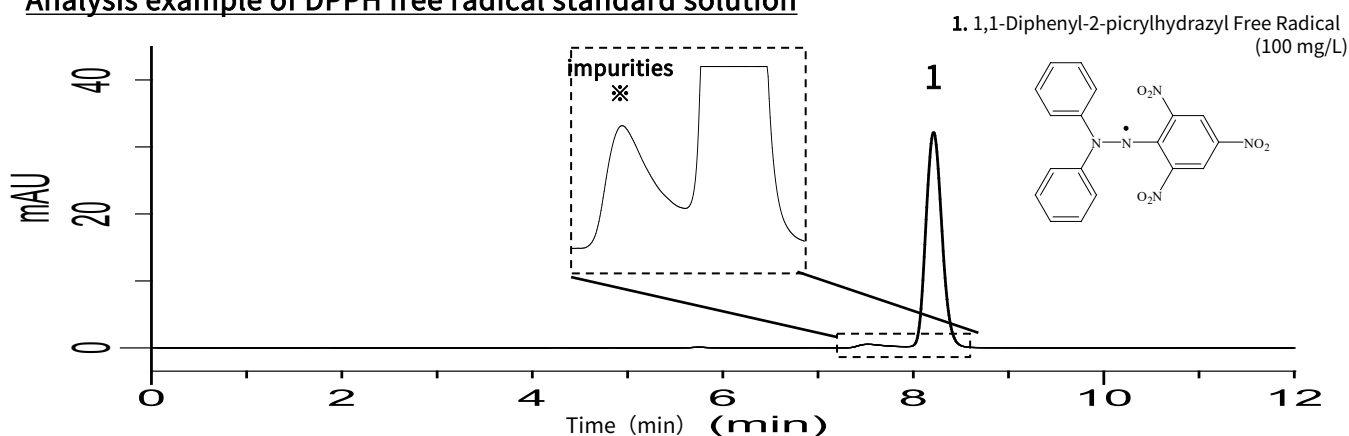


Analysis example of DPPH standard solution

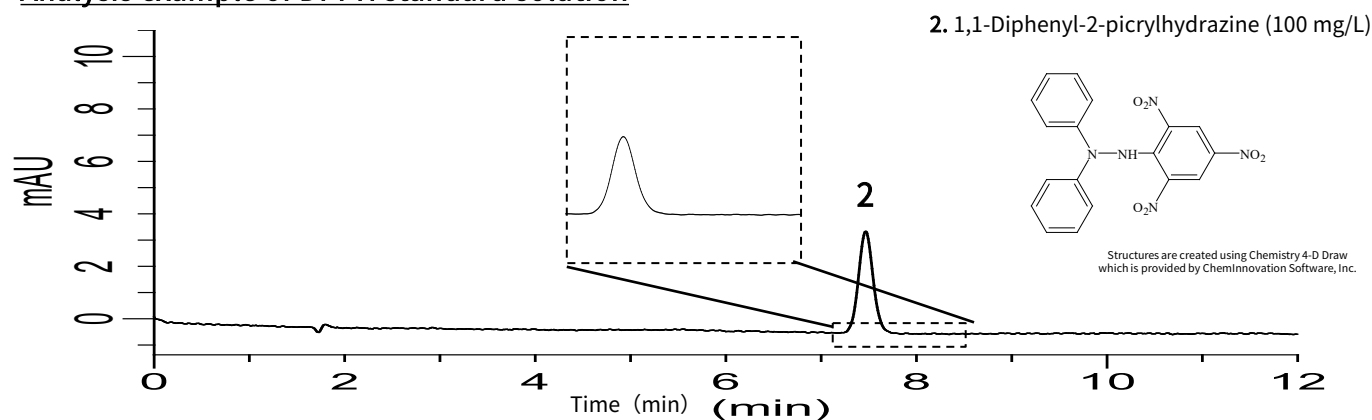



HPLC Condition

Column : InertSustainSwift C8 (5 μ m, 150 $\times$ 4.6 mm I.D.)
 Cat. No. : 5020-88331
 Eluent : A) CH₃CN
 B) H₂O
 A/B = 55/45, v/v
 Flow rate : 1.0 mL/min
 Col.Temp. : 40 $^{\circ}$ C
 Detection : VIS 540 nm (DAD)
 Inj.Vol. : 10 μ L

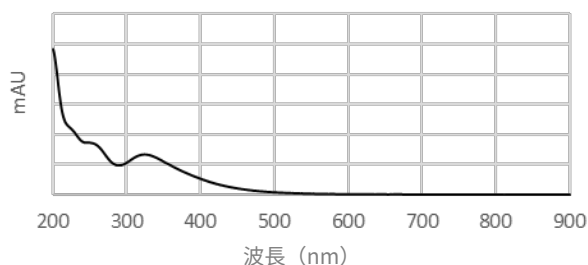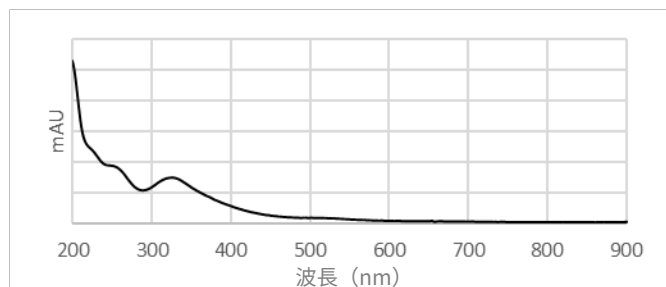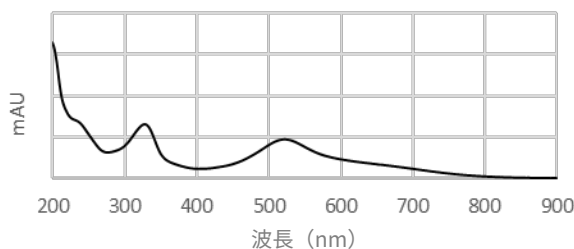

An example of analysis of DPPH free radicals at low concentrations is shown. Impurity (*) peaks appear to be tailing, but since no tailing was observed with the standard product on the previous page, it is presumed to be decomposition within the HPLC.

Analysis example of DPPH free radical standard 10 mg/L

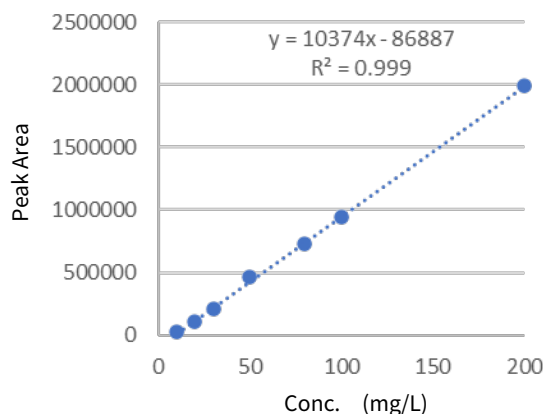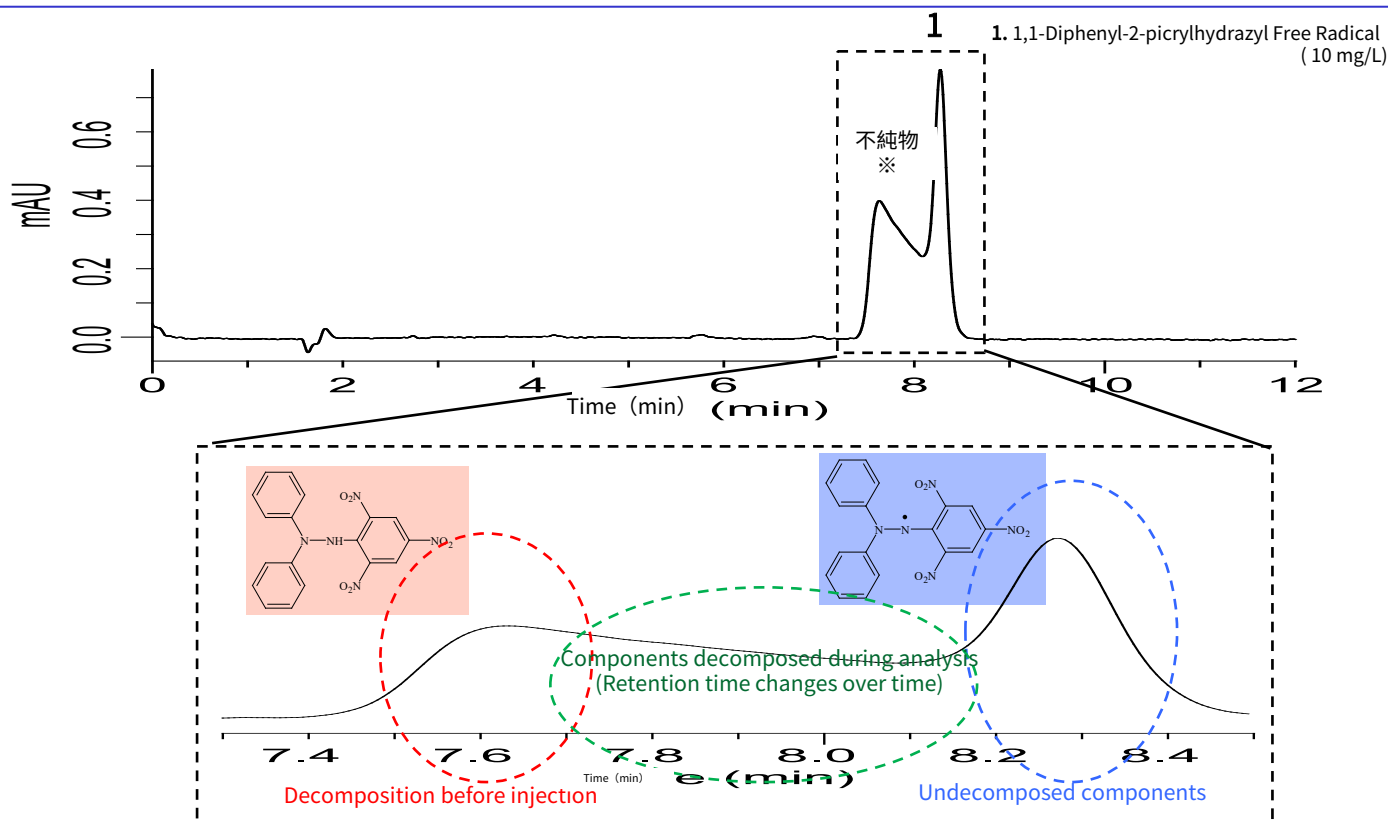

DPPH free radical calibration curve

GL Sciences disclaims any and all responsibility for any injury or damage which may be caused by this data directly or indirectly. We reserve the right to amend this information or data at any time and without any prior announcement.

GL Sciences Inc. Japan

22-1 Nishishinjuku 6-chome
Shinjuku-ku, Tokyo
163-1130, Japan

Phone: +81-3-5323-6620
Fax: +81-3-5323-6621
Email: world@glsc.co.jp
Web: www.glsciences.com

GL Sciences Inc. USA

4733 Torrance Blvd. Suite 255
Torrance, CA 90503
USA

Phone: +1-310-265-4424
Fax: +1-310-265-4425
Email: info@glsciencesinc.com
Web: www.glsciencesinc.com

GL Sciences B.V.

Dillenburgstraat 7C
5652AM, Eindhoven
The Netherlands

Phone: +31-40-254-9531
Email: info@glsciences.eu
Web: www.glsciences.eu

GL Sciences (Shanghai) Limited

Tower B, Room 2003
Far East International Plaza
No.317 Xianxia Road, Changning District
Shanghai, China 200051

Phone: +86-21-62782272
Email: contact@glsciences.com.cn
Web: www.glsciences.com.cn



International Distributors

Visit our Website at www.glsciences.com/distributors